


Quality guaranteed
Outstanding performance



injection



surgery

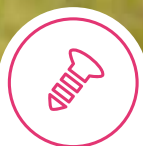

orthopedics



hospitalisation



capture



protection



hygiene




livestock



grooming

génia AUSTRALIA



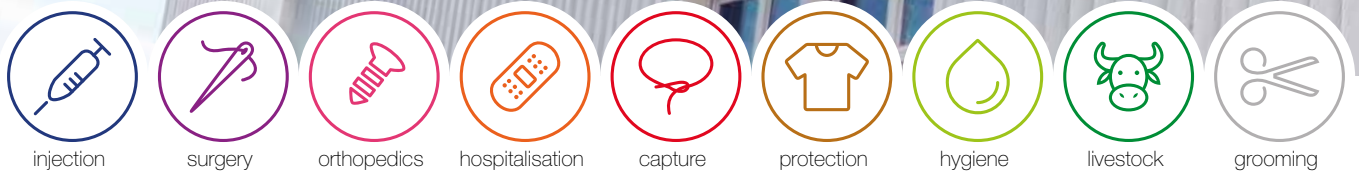

INTRODUCING GENIA FRANCE AND GENIA AUSTRALIA

Génia

Since its founding in 1945, Génia has been at the forefront of innovation, driving its growth and establishing itself as a market leader in providing specialised solutions to veterinarians worldwide.

Headquartered in France, Génia is a family-owned company with over seventy years of expertise in the veterinary supply industry.

Our journey began with the development of the pioneering Plexiglass syringe called Géniaplex, and today, we manufacture +4000 SKUs across nine product ranges, available in over seventy countries. With branches in the UK, Australia, and the USA, we are able to stay ahead of global market trends, continuously delivering high-quality products that elevate veterinary practices with a focus on animal comfort and hygiene.



Génia Australia

For more than 13 years, Genia Australia has been an integral part of the Australian veterinary landscape, delivering innovative solutions and high-quality products to veterinary professionals across the country.

Since our establishment in 2011, we have worked tirelessly to understand and meet the unique needs of Australian veterinarians and nurses, building lasting relationships with clinics, wholesalers, veterinary groups, and industry partners.

We always make it a priority to visit veterinary clinics regularly, ensuring that we maintain a direct connection with the professionals who use our products every day. This hands-on approach allows us to listen to their needs, gain valuable insights, and continuously improve our offerings.

We work closely with the leading wholesalers and some veterinary groups to ensure that Genia products are easily accessible to those who rely on them.

We also actively participate in conferences, where we showcase our products and engage with the wider veterinary community. These events provide us with an opportunity to share our latest innovations as well as a good “baguette” !

Our new motto, **«The Frenchy Team aligning with Aussie Legends»** reflects our strong commitment to working in harmony with the incredible vets and nurses who make up the backbone of the Australian veterinary industry.

It's always a privilege to collaborate with such dedicated professionals who share our vision of improving the comfort and well-being of animals. We believe that working with passionate, skilled individuals is what truly drives our success and enables us to continue delivering solutions that enhance veterinary practices.

At Genia Australia, we are proud to be part of a global family dedicated to advanced animal healthcare, and we look forward to the many opportunities ahead to support the vets and nurses who care for our beloved animals.

And if ever you would like to speak in French or simply say «bonjour,» we are here for you!

Contact us:

www.geniaaustralia.com

contact@geniaaustralia.com



Hospitalisation	p.8
Surgery	p.22
Orthopedics	p.25
Protection	p.26
Injection	p.28



Surgery	p.38
---------------	------



Warming system	p.44
Blood pressure	p.47
Electrosurgery	p.50



Enrichment toys	p.56
Feeding	p.59
Rest	p.62
Apparel	p.63
Walking accessories	p.65
Hygiene	p.66



Protection	p.70
Injection	p.72
Surgery	p.74
Livestock	p.75
Hospitalisation	p.80



SMALL ANIMALS

SURGICAL INSTRUMENTS

EQUIPMENT

vétopop
The best for my pets

LARGE ANIMALS & EQUINE




AUSTRALIA
génia

KEEP CALM AND GO GREEN

**BE PART OF THE
SOLUTION**

with Génia sustainability commitment

VELFAST

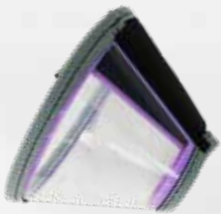

FUNCOL



TRANSCOL



GREENCOL



GREENCOL
MADE OF
80%
RECYCLED
MATERIAL

NEOCOL



ALL E-COLLARS
100%
RECYCLABLE



RECYCLING

ELIZABETHAN COLLARS

P. 8-11



TIKER

P.19

TIKER®

100%
RECYCLABLE



**I'm
green™**
Plastic
98% Biobased
Verified by ASTM D6866



RECYCLING



**BIODEGRADABLE
SLOW FEEDER BOWLS**


P. 44-45



ORGANICS

 **KEEP CALM AND** 



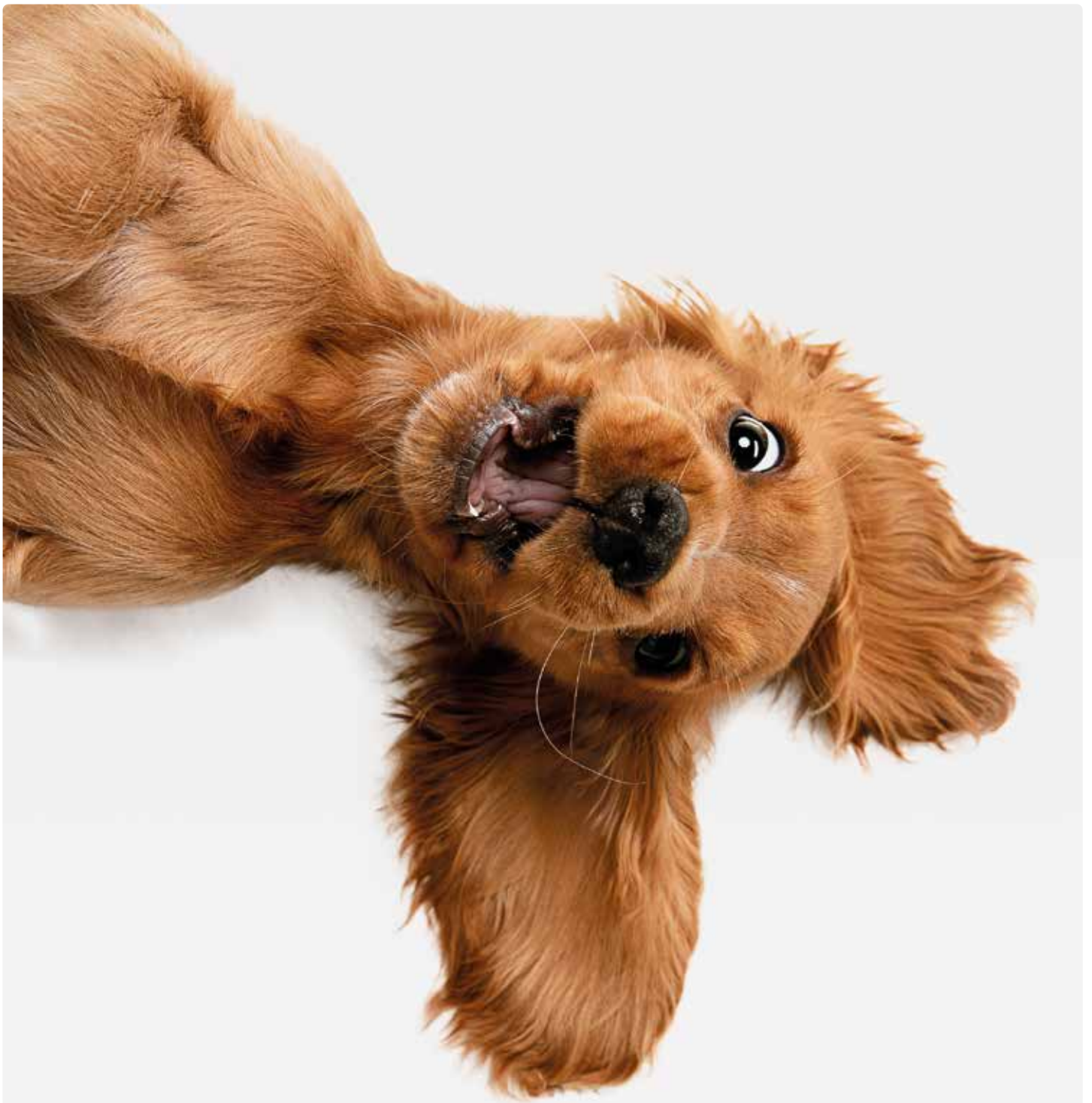

et voila !



SMALL ANIMALS



Our range of Elizabethan Collars has been carefully thought out to fit all vet practices needs, thanks to the following features:

Our well-known quick fastening system: Release the E-Collar by bending the inside part completely back, until it marks white.

Scalloped edge: The wavy outer edge makes them sturdier and less sharp on your legs.

Smart hanging system: The open storage hole allows you to store each E-Collar on the same hanger; pull off just the one you need. Save space and time!

Colour coded according to the size: For better recognition, making your life easier.

3 to 5 neck adjustments: To get the perfect snug fit for any animal.

ALL our E-Collars are 100% recyclable.

Code Transcol	Code Funcol®	Code Greencol	Code Neocol®	Code Velfast®	Size
306100	306200	306300	306420	306450	7.5 cm
306102	306202	306302	306422	306452	10 cm
306104	306204	306304	306424	306454	12.5 cm
306106	306206	306306	306426	306456	15 cm
-	306218	306307	-	-	17.5 cm
306108	306208	306308	306428	306458	20 cm
306110	306210	306310	306430	306460	25 cm
306112	306212	306312	306432	306462	30 cm
-	306213	306314	-	-	35 cm
-	306214	306315	-	-	40 cm

Recyclable E-Collars

VELFAST®



GREENCOL®



FUNCOL®



TRANSCOL®



NEOCOL®



ECO-FRIENDLY

GREENCOL

GREENCOL is our economical and ecofriendly E-Collar.

Made from 80% recycled transparent polypropylene and 100% recyclable.

Ideal for clinics that are trying to reduce their carbon footprint by using a cost effective range of E-Collars.

An "IN-BETWEEN" sizes is now available: 17.5 cm

- Transparent polypropylene: 80% recycled, 100% recyclable.
- Quick fastening system.
- 3 to 5 neck size adjustments.
- Wavy outer edge.

Sold individually or in boxes of 12 e-collars.



NEW SIZE
17.5 CM

80 % RECYCLED
100 % RECYCLABLE



VELFAST®

VELFAST® is our quickest and easiest E-Collar. With its strong Velcro fastening system, VELFAST® is a time saver! It also features soft double edges, for more comfort to the animals neck.

- Transparent polypropylene, with coloured neoprene edges. Recyclable.
- Colour-coded according to size.
- Velcro quick fastening system.

Sold individually or in boxes of 12 E-Collars.



QUICK &
EASY



FUNCOL®

FUNCOL® is exactly what it says: it is fun!

Life is better with **COLOURS**. The coloured paw prints will make your team and customers smile.

- Coloured paw prints on transparent polypropylene, 100% recyclable.
- Colour-coded according to size.
- Quick fastening system.
- 3 to 5 neck size adjustments.
- Wavy outer edge.



FUN I

An in-between sizes is available: 17.5 cm

If the 15 cm is too small and the 20 cm is too big you can now use the 17.5 cm for a perfect fit!

Sold individually or in boxes of 12 e-collars.

ADDITIONAL
SIZES

17.5 CM

35 CM

40 CM




TRANSCOL

TRANSCOL is our 'chic' E-Collar.

It gives a French touch to your four legged friends!

- Translucid coloured polypropylene, 100 % recyclable.
- Colour-coded according to size.
- Quick fastening system.
- 3 to 5 neck size adjustments.
- Wavy outer edge.

Sold individually or in boxes of 12 e-collars.

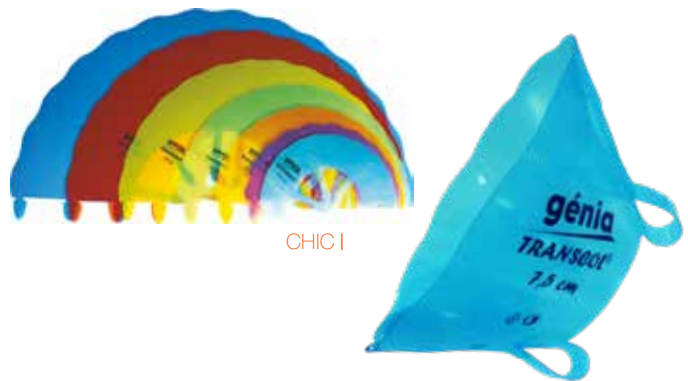

CHIC I

NEOCOL®

NEOCOL® has a comfortable shock absorbing neoprene outside edge.

Specially made to provide comfort to both patients and owners. Softer on owners legs and protects furniture & walls.

Feel safe on your night walks, thanks to its reflective layer.

- Transparent polypropylene, with a coloured neoprene edge, recyclable.
- Colour-coded edge according to size.
- Quick fastening system.
- 3 to 5 neck size adjustments.

Sold individually or in boxes of 12 e-collars.



COMFORTABLE I

DOG E-COLLARS HANGER

All our E-Collars have a patented hanging system. Store them all on one hook and just pull off the one you need. It saves you both time and space.

Code Product

306117 Hanger



DESIGNED TO FIT
ALL GENIA E-COLLARS



FUNCOL®

MANJIMUP VET CLINIC, WA

Shannon, Vet Nurse

Genia E-collars are perfect, easy to put together compared to others brands.

RINGWOOD VET CLINIC, VIC

Briar, Practice Manager

We did get the 17.5cm e-collars in, and they have been working a treat.

SWEET PEA ANIMAL HOSPITAL, NSW

Samantha, Practice Manager

We switched to Funcol E-collars because we like the paw prints pattern, we found them nicer than plain ones. We never had any issue with them so far!

TESTIMONIALS

GREENCOL

HAMILTON HILL VET HOSPITAL, WA

Bruce, Vet Nurse

The Greencol E-collars are fantastic. They are pretty easy to put together and are durable.

ULVERSTONE VET CLINIC, TAS

Shanna, Vet Nurse

The Greencol E-collars are cheaper than Buster. We also like the fact we can hang them on the same hook, and pull the one we need without taking off the others ones, thanks to the opened hole.

VETS ON BALWYN, VIC

Mikaela, Practice Manager

We like the Greencol E-collars because they are easy to put together and release. It's ideal when you need to resize them on the patient.

BRIGHTON VETERINARY SERVICES, TAS

Karinda, Vet Nurse

We find the Greencol E-collars more comfortable than everything, they are softer on the neck. They are easy to put too.

TESTIMONIALS

NEOCOL®

EASTERN SHORE, TAS

Kea, Practice Manager

I like Neocol E-collars because of their soft edges and the fact they are reflective. They are also easy to pull apart when they are in the wrong side.

LARA PET HOSPITAL & VETERINARY CENTRE, VIC

Christie, Vet Nurse

All the staff love the Neocol E-collars! We like the soft padded edge and find these e-collars more flexible and comfortable than others brands.

ARMSTRONG CREEK VET PRACTICE, VIC

Rachel, Vet Nurse

We love the Neocol E-collars, they are easy to click! The soft edge protects the patient from scratching himself.

KUDAH VETERINARY SERVICES, WA

Amy, Vet Nurse

We like the Neocol E-collars because they are user-friendly; they're a nicer for your legs when pets run into you.

DANDENONG RANGES VET CLINIC, VIC

Elissa, Vet Nurse

We use the Neocol mostly for smaller animals (such as puppies or kittens) or for scared patients, as it is more comfortable for them.

TESTIMONIALS

TESTIMONIALS

TRANSCOL

THE ANIMAL HOSPITAL KALGOORLIE, WA

Marge, Practice Manager

We like the colourful Genia E-collars. They are quite sturdy and last really well.

VELFAST®

CLAREMONT VETERINARY SURGERY, TAS

Kathryn, Vet Nurse

The Velfast E-collars are really convenient. They get on and off quite easily. I also like the soft edge on the neckline. The Velcro is strong, making the E-collar re-usable, compared to a clickable E-collar, less sturdy. They look nice too.

PAUL'S VETERINARY CLINIC, TAS

Jess, Vet Nurse

The Velfast E-collars are easy to put on. They're ideal for patients with small neck, as we can easily adjust them; not too loose or too tight. They're comfortable for the animals. Nice and easy.

CANNING VALE VET, WA

Sophie, Vet Nurse

Just got some more, we absolutely love these e collars. They are so easy to use and even clients comment on the look and durability they have! Thanks so much for recommending them to us.

ANIMAL WISE TRAINING, VIC

Dr Caroline, Vet

I have done some initial trials, in a home setting, to help my little old miniature pinscher, Roo, get used to the Velfast E-Collar, and I was very impressed how quickly he came to accept it! Like most dogs, he has not been a fan of E-Collars. It's amazing what taking your

CASEY'S BEACH VET CLINIC, NSW

Christine, Head Vet Nurse

We love the VELFAST, especially on angry kitty cat, and our customer do too.

MELBOURNE EYE VET, VIC

Tracey, Vet Ophthalmic Nurse

Dogs tolerate more the Velfast E-collars than any others e-collars, as they are more comfortable. Owners are surprised by this product, and grateful because their pets are less stressed.

CAPE HAWKE VET HOSPITAL, NSW

Cass, Nurse

These VELFAST, are so easy to use. We had one patient who was able to remove the normal e-collar with his feet. He has unable to do this to the VELFAST.

MIDWAY ROAD VET CLINIC - ELIZABETH EAST, SA

Georgia, Head Nurse

Genia Velfast collars are our favourite collars to use post-operatively. The padding on the collars are comfortable for patients as well as owners and the velcro makes them easy to place on all pets.

SALISBURY VET CLINIC, SA

Bec, Head Nurse

The GENIA Velfast collars are amazing, a game changer, time saving product for the vet industry. So much more comfortable for the patients too.

CROOKWELL VET HOSPITAL, NSW

Dr Jayde, Senior Veterinarian and Director

They are easy to put on, there is no mocking around with the Velfast collar and they stay on

TESTIMONIALS



VELFAST®



CAROLINE SPRINGS VET - Dora and Peanut - VELFAST®



ALL ANIMAL EYE SERVICES - Jet - VELFAST®



NORTHSIDE VETERINARY CENTRE
Kristofferson - GREENCOL



ADVANCED VETCARE - Maggie and Merlin - GREENCOL



GREAT WESTERN ANIMAL HOSPITAL - Evie - GREENCOL



LENNOX HEAD VET CLINIC - Lizzie - FUNCOL®



SUPERVETS - Kate The Vet - FUNCOL®



WILLOUGHBY VET HOSPITAL - Sunny - FUNCOL®



VETWEST - Yoda - TRANSCOL



NORTH COAST VETERINARY SPECIALIST
AND REFERRAL CENTRE - TRANSCOL



SASH VETS - Andy The Adventurer - TRANSCOL



Oogie the Ruby Cavoodle - NEOCOL®



STEPHEN TERRACE VET CLINIC - Poppy - NEOCOL®



RINGWOOD EAST VET CARE - Bax - NEOCOL®

STRAP BAND

Cohesive elastic bandages.

Available in 5 cm, 7,5 cm and 10 cm.

- Economical self-adhesive bandages.
- Double layered: strong and sweat / water-resistant.
- Easy to tear: no scissors needed, rip in a straight line.
- Coloured starting strip: find the edge easily.
- Purple bitter taste cohesive bandages.

Available in 5 cm, 7,5 cm and 10 cm.

- Strap band with a bitter taste to prevent from chewing, biting and tearing off the bandage.
- Contains denatonium benzoate.

Singly wrapped. Sold individually or in boxes of 24 rolls for size 5 cm and boxes of 12 rolls for the 7,5 cm and 10 cm.

Blue	Red	Green	Yellow	Black	Bitter taste	Size	Pack.
307202	307200	307201	307203	-	307230	5 cm x 4.5 m	x 24
307208	307206	307207	307209	-	307232	7.5 cm x 4.5 m	x 12
307214	307212	307213	307215	307216	307234	10 cm x 4.5 m	x 12

Code	Content
307220	Kit of 32 Strap Band: 5 cm x 4.5 m : 2 x 4 x 7.5 cm x 4.5 m : 2 x 4 x 10 cm x 4.5 m : 2 x 4 x
307222	Kit of 16 Strap Band Bitter taste: 5 cm x 4.5 m 8 x 7.5 cm x 4.5 m 8 x



BITTER TASTE BANDAGES

VICTORIA PARK VET CLINIC, WA

Monique, Head Vet Nurse

Bitter taste bandage gives us an extra security and more comfortable. It ensures the patient doesn't lick it and secures the catheter and the fluid.

WALLAN VET HOSPITAL, VIC

Janelle, Practice Manager

Genia bitter taste bandage is the best bitter one we've tried; it's successful. We use it over the catheter.

TESTIMONIALS

STRAP BAND

BUNBURY VET CLINIC, WA

Amanda, Head Vet Nurse

Genia bandages are easy to rip. We like as well the little marker to find the edge; it's a big help!

BANGALOW VETS, NSW

Flo, Head Vet Nurse

The quality of Genia bandages is really good for the price.

FUN STRAP®

PETS ♥ THEIR VETS. As they are not able to say it, we wrote it down for them on our FUN STRAP®. Available in 5 cm, 7,5 cm and 10 cm.

- Funny self-adhesive bandages.
- Latex-free.
- Strong, sweat and water-resistant.
- Easy to tear: no scissors needed, rip in a straight line.
- Coloured starting strip: find the edge easily.

Singly wrapped. Sold in boxes of 24 rolls for size 5 cm and boxes of 12 rolls for the 7,5 cm and 10 cm.



NEW
DESIGN

LATEX
FREE

MIXED
PACK !

Code	Size	Content	Pack.
307240	5 cm x 4.5 m	8 rolls of each design	x 24
307241	7.5 cm x 4.5 m	4 rolls of each design	x 12
307242	10 cm x 4.5 m	4 rolls of each design	x 12
307243	Set - 3 rolls of each size and each design		x 27



SOFT SPLINT ROLL

These splints can be moulded to the shape you need.

Lightweight, compact and very comfortable for the animal: avoid pressure sores.

- Reusable: easily reshapable for another application.
- Washable: closed-cell foam allows easy cleaning or disinfection (material compatible with most cleaning products).
- Easy to cut.
- Radiolucent.
- Soft and resilient: recommended for temporary immobilisation.
- Ideal for any small animals as well as foals and goats.
- 3 sizes, colour coded.

Singly wrapped. Sold individually.



SHAPEABLE & RADIOLUCENT

Code	Dimensions
3078301	4.5 x 25 cm
3078321	5.5 x 35 cm
3078341	9 x 40 cm

GENIA GUTTER SPLINTS

Plastic gutter splints ideal for fractures immobilisation. Comfortable, non-irritating thick foam padding and multiple ventilation holes.

- Very strong material, poly-gloss finish.
- No sharp angle, thick closed-cell protective foam that extends beyond the edges of the splint to protect the animal.
- Can be cut to adjust to the animal's leg size.
- Air holes along the length of the splint allow air to circulate.

Singly wrapped. Sold individually.

LIGHT WEIGHT & BREATHABLE



Code	Length
307801	20 cm
307806	30 cm
307811	40 cm

TESTIMONIALS

GUTTER SPLINTS

DARDANUP VET CLINIC, WA

Dr Liz, Vet

Genia gutter splints have a good shape and fits perfectly the patients, especially dogs. We like the foam padding.

FLAT SPLINTS

Splints with an aluminium side and a protective soft foam side. Ideal for tails.

- Can be cut.
- Radiotransparent.
- Sturdy, light and atraumatic.
- Ideal for a temporary immobilisation.
- Available in 4 widths.

Singly wrapped. Sold individually.



PRICE DROP

Code	Length
3078201	2 x 20 cm
3078241	5 x 40 cm
3078261	8 x 45 cm
3078281	10 x 45 cm

TESTIMONIALS

SOFT SPLINT ROLL

LEEMING VET CLINIC, WA

Amy, Vet Nurse

We used the soft splint rolls mainly on small dogs, such as pugs; it went really well! We can cut the splints into the length required, which is very convenient.

WEST QUEANBEYAN VET HOSPITAL, NSW

Katherine, Vet Nurse

We do lot of orthopaedics and we found that the cat tolerate the soft splint a lot better.

HERITAGE VET CLINIC, VIC

Dr Rob, Vet

They are excellent for particular cases, especially for small dogs or cats. They're light and work really well.

MARIBYRNONG VET CLINIC, VIC

Jacqueline, Practice Manager

I like the fact we can cut the soft splint rolls to shape and mold them.



COSYPAD®

Non-absorbent meshed mat providing comfort and hygiene; it has become a clinic favourite!

💧 The Cosypad® is non-absorbent, allowing liquids to go through the mat, leaving the animal dry. Perfect to use in your cages for incontinent patients or when waking up from GA's. It minimises pressure on muscles, joints, and prevents bedsores, making it ideal for hospitalised patients.

🔥 Thanks to its anti-slip layer, the Cosypad® is perfect on your consult or surgery tables. The meshed material provides optimal comfort and ventilation, as well as a thermal barrier to cold tables.

🦷 The small Cosypad® makes a great dental mat: place the head of the animal on top of the mat while descaling, and let liquids drain through.

😴 Tired ? Use your Cosypad® as an anti-fatigue mat !

Last but not least...

☀️ Easy to clean: Cosypad® is machine washable (up to 40°C) and dries quickly even in winter.

🔄 Chew & stain resistant: Cosypad® is designed to last years !

IDEAL FOR YOUR:

Consultation
tables

Surgery
tables


Dental
procedures

Cages

Code	Size	Use
327200	33 x 50 cm	Ideal for dental
327201	39 x 71 cm	Ideal for cages / consultation tables
327203	55 x 71 cm	Ideal for cages / consultation tables
327205	60 x 115 cm	Ideal for bottom cages / floor

♻️ Recyclable mat

HANGING
HOLE

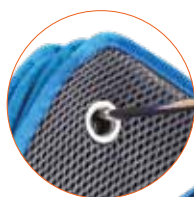

ANTI-SLIP LAYER
UNDERNEATH



MACHINE
WASHABLE !



40°C - Do not tumble dry
Wash in delicate cycle
if using top loaders

ROUNDED
CORNERS



ROSEWOOD VET SERVICES, QLD

Jodi, Head Nurse

We use the Cosypads for cages, they provide padding and comfort. Easy to clean and non-slip. Would 100% recommend.

SYDNEY ROAD VET CLINIC, NSW

Dr Kayleigh, Vet

Our Practice Supervisor Piper adores her pink Cosypop. She likes to use it as a bed and playpen, and so she can be transported within areas of the clinic with ease. We love that it is almost our clinic colours of magenta and grey. Kay Gerry BVSc, Sydney Road Veterinary Clinic, Balgowlah.

TESTIMONIALS



DENISTONE VETERINARY SURGERY - Tong Tong - COSYPAD®

Perfect to protect from the cold of consultation tables

SUPERVETS, QLD

Dr Anthony, Vet

"As a veterinary nurse for 12 years before a practicing vet, I can thoroughly appreciate the sheer amount of cleaning, washing and disinfecting required in day to day operations of a hospital. We want our patients to be comfortable and warm and are always looking to improve hospital care and efficiencies when dealing with bodily excrements.

When Génia introduced their new COSYPAD® we trialed them, and within weeks filled every cage in 2 24hr hospitals with them as routine. Now every patient has a comfortable non slip bed which keeps them off cold surfaces and prevents urine scalding. Being hygienically washable means we have drastically reduced our use of towels and the costs of washing as a consequence. We have also found them particularly useful as a non-slip surface in consult and prep benches, and routinely use them as the resting surface for patient heads during dental procedures, with pee pads underneath keeping their face far cleaner. Very useful in the car too ! We absolutely love them! :)"

TESTIMONIALS



DENISTONE VETERINARY SURGERY - Bella and Dougal - COSYPAD®



Ideal for dental surgeries

MCINTYRE VET SERVICES, VIC

Sally, Practice Manager

We use our Cosypads everyday, either for dentals or on our surgery table. They're really good, we're happy with them!

THE CAT CLINIC HOBART, TAS

Michelle, Head Vet Nurse

We use the small Cosypad mostly for our dental procedures. As the liquids go through, it keeps the cat dry. We like to use it as well on our consultation tables, during vaccinations or when taking blood samples. The mat doesn't move thanks to its anti-slip layer.

KARRINYUP SMALL ANIMAL HOSPITAL, WA

Emma, Head Vet Nurse

I'm really liking the cosypads for our dental procedures. They provide much more comfort for our patients than what we were previously using.

DROUIN VET CLINIC, VIC

Linda, Practice Manager

We like that the Cosypad is machine washable. The soft padding makes it ideal for orthopaedic recoveries.

LYNWOOD VET CLINIC, WA

Dr Louise, Vet

The Cosypad padding is really good as it avoids pressure sores. We use the mat also for our dental procedures.

CANNINGTON VET, WA

Dr Trevor, Vet

We value the Génia Cosypad for its efficiency in keep our patients on a dry comfortable surface but also for its hardwearing quality and ease of cleaning. Would recommend as a good investment.

WAUCHOPE VETS, NSW

Sarah, Vet Nurse

We use our Cosypad for all our dentals. As the liquids go through, the patient doesn't sit in a cold wet surface. And it's easy to wash! We like as well its non-slip layers.

COHUNA BORDER VETS, VIC

Kelly, Vet Nurse

We use the Cosypad in our cages. As the liquids go through, the patients don't stay in their own wee. We use it as well on our examination table because it has anti-slip layers.

NORTH COAST VETERINARY SERVICES, NSW

Dr Lauren, Vet

We are currently using the Cosypad for a dog with snake bite paralysis, and they are going really well.

TAREE VET HOSPITAL, NSW

April, Practice Manager

The Cosypads are great, we have been using them for dentals, also older dogs and longer staying patients in hospital.

TESTIMONIALS



VACUUMAT®

VACUUMAT® - Our lightweight positioning mat that holds patients in the desired position prior to a surgical procedure, ultrasound, MRI, CT or X-ray (as it is radiolucent). Its insulation and shaping properties are achieved through microbeads made of sturdy polyurethane material.

You can easily clean it with soapy water or simply wipe it down.

The material is thick and resistant to damage.

It comes as a complete kit, including the mat and an adapted vacuum pump.

How does it work ?

- 1. Lay the animal on the mat, in the desired position.
- 2. Close the valve.
- 3. Suck the air out of the mat with the adapted pump in less than 10 seconds.

Et voilà!

Why do you love it?

- The mat conforms to the body shape and holds the patient in place. The animal won't move or slide away.
- It is ideal for X-rays positioning, especially for hips and legs, that often require difficult angles.
- As it provides stable immobilisation, the VACUUMAT® is also perfect for eye surgeries.
- The mat gets stiff while remaining flexible and very comfortable for the animal. The stiffness can be adjusted by slowly releasing the valve.

What else?

- The VACUUMAT® reduces hypothermia as it keeps patients off cold surfaces, and it is highly insulating. The microbeads hold and maintain patients' body heat, while preventing from a drop in temperature.

THERMAL
INSULATION

EASY TO CLEAN

POSITIONING
MAT



Code	Description	Size
327302	VACUUMAT® with pump included	40 x 60 cm
327304	VACUUMAT® with pump included	60 x 100 cm
327306	Additional pump	-

VACUUMAT

MIDWAY ROAD VET CLINIC - ELIZABETH EAST, SA
Georgia, Head Nurse

We love using our Genia Vacuumat during anaesthetics. It is a great positioning aid which also helps the patients hold a regular temperature and is super easy to clean!

MELTON VET CLINIC, VIC
Melissa, Practice Manager

We love using the Genia Vacuumat at our clinic. It moulds around our patients keeping them in a secure position during procedures. It is easy to use and clean, and helps to keep our patients warm. We are very happy with the product.

TESTIMONIALS



**TASMANIAN ANIMAL HOSPITALS - BELLERIVE, TAS***Claire, Practice Manager*

I love the Vacuumat to position the patient on the surgery table. I use it every day.

HUON VALLEY VET HOSPITAL, TAS*Sarah, Vet Nurse*

The Vacuumat is quite good. We use it mostly for X-rays procedures and for tall animals as it holds their position.

LINDISFARNE VET CLINIC, TAS*Dr Sam, Vet*

We find the Vacuumat convenient to mold around the patient body. We use it mostly to stabilise the animal on his back during surgery.

ROUSE HILL FAMILY VETS, NSW*Belinda, Practice Owner*

We use the Vacuumat to assist in comfortable positioning for patients during ultrasounds and radiography procedures.

GREENCROSS VETS EARVILLE, QLD*Emily, Practice Manager*

We love the vacuumat, even with a little bit of blood it comes off very easily. The pump works fine. It's also such a quick-release that you can shape it to your needs. It really molds around the patients limbs and you can adjust the firmness. We love it for dentals, X-rays, surgeries, we basically use them all the time!

BICTON VET CLINIC, VIC*Dr Kay, Vet*

We use the Vacuumat for eye surgeries, as it maintains the patient on the correct orientation all along the operation.

WELLINGTON VILLAGE VET, VIC*Kelly, Vet Nurse*

The Vacuumat is easier than sand bags. It's ideal to use for the vets, to position the patient.

DEEPWOOD VET HOSPITAL, WA*Michelle, Vet Nurse*

The Vacuumat is soft for the animals and keeps them in the position that we want.

TESTIMONIALS

CANNINGTON VET, WA*Dr Trevor, Vet*

The Genia Vacuumat is so helpful and versatile for patient cushioning, insulation and positioning. We have been using them every day for surgery and radiographs for a couple of years now. It is easy to clean and showing no sign of wear. I can confidently recommend them.

EXOTIC VET NURSE TRAINING, WA*Samantha, Owner*

One challenge we face that I think we can all agree on keeping patients in position during procedures, especially the deep chested dogs! I've been looking into the best options for this, as towels, tie downs and foam pads just weren't cutting it anymore, but I've since come across the Vacuumat from Genia Australia that might be a one size fits all.

The mat is insulated to maintain body heat and inflates in less than 10 seconds, around your patient in their desired body position. No more rolling and folding of towels and repositioning half way through surgery when they slip.

THE DOG CLINIC, TAS*Dr Raj, Vet*

We love the the Vacuumat, we have one small and one large that we use mostly for our desexing procedures.

ANIMAL MEDICAL CENTRE, TAS*Alice, Head Vet Nurse*

We use our Small and Large Vacuumat everyday for all our surgeries. On top of positioning the patient, they transfer the heat, maintaining his body warmth. We grab them for X-rays as well. We simply love them!

**YARRAVILLE VILLAGE ANIMAL HOSPITAL, VIC***Steph, Vet Nurse*

We find the vacuumat helpful during surgery in keeping animals positioned, especially for any midline surgical procedures.

LARKHILL VETS, WA*Hayley, Head Vet Nurse*

The Vacuumat makes the animal more comfortable. It's way better to position him. It's simply awesome! We use it for both surgery and X-rays.

CAMBERWELL VET CLINIC, VIC*Vivian, Vet Nurse*

The Vacuumat is a big help, especially for deep chested dogs such as greyhounds! We mostly use it for our CT scans.

VICTORIA PARK VET CLINIC, WA*Monique, Head Vet Nurse*

We have the Vacuumat in 2 sizes; it's much better than using sand bags!

PORT MACQUARIE VET HOSPITAL, NSW*Dawn, Practice Manager*

We like how the Vacuumat is easy to use and to store. It's a reliable product.

BRIMBANK VET CLINIC, VIC*Merinda, Practice Coordinator*

The Vacuumat is quite handy; we use it mostly for our X-rays as we don't need to be in the room anymore to hold the patient. We use it as well for few surgeries, especially on rabbits.

TESTIMONIALS





HOSPITALISATION

DIGIFLASH

Digital thermometer.

- Flexible tip.
- Completely waterproof.
- Records the last measurement.
- Shock resistant.
- Temperature measurement in 10 seconds.
 - Temperature range: 32 to 42 °C .
 - Accuracy: 0,1 °C.
 - Alarm indicating measurement end.
- Low battery symbol display.
- Auto power-off.

Sold individually in a transparent polycarbonate box.

Code	Product	Tip length	Total length
227972	Digital thermometer	75 mm	155 mm



QUICK READING: |
10 SEC ! |



DIGIFLASH

SOUTH PAWS MORNINGTON, VIC

Jennine, Vet Nurse

The Digiflash thermometers are quick and we like the flexible tip. It's also water proof, I left it in the washing machine and it kept working!

MOUNT EVELYN VET CLINIC & ANIMAL HOSPITAL, VIC

Kellie, Vet Nurse

We find the Digiflash way easier to use thanks to its flexible tip. It's also more comfortable for the patient.

TESTIMONIALS

DIGIFIX

Digital thermometer.

- Completely waterproof.
- Temperature measurement in 10 seconds.
- Records the last measurement.
 - Temperature range: 32 to 43°C.
 - Accuracy 0,1°C.
 - Alarm indicating measurement end.
- Thinner, ideal for exotics
- Low battery symbol display.
- Auto power-off.
- Shock resistant.

Sold individually in a transparent polycarbonate box.

Code	Product	Tip length	Total length
227978	Digital thermometer	75 mm	155 mm



10 SEC |

NURSE WATCH

Tick Tock! It's time to keep calm and love your nurse!

- Waterproof silicone watch.

Code	Model
467600	Nurse watch



TOURNIQUET GARFLEX

Stretch-fabric quick-release tourniquet.

Sold individually.



Code	Size
307860	390 x 25 mm



TIKER®

Reusable, quick and easy to use, our tick remover is an essential tool, for safely removing small and large ticks. Our special spoon design helps you to remove ticks without leaving mouthparts in the skin or compressing the tick's abdomen. Minimises the transfer of infectious agents.

- Made from Green Polyethylene, the TIKER® is eco-friendly: the bio-based material makes it 100% recyclable.
- The cardboard packaging is 100% recycled and recyclable.

Sold individually or in boxes of 24 tick removers.

Code **Model**

191061 Tiker®



I'm green™
Plastic
98% Biobased
Verified by ASTM D6866



ECO-FRIENDLY
BIO-BASED MATERIAL
100 % RECYCLABLE

DRYSOFT®

Economical and highly absorbent under pads.

Made of a SAP polymer gel layer (Sodium polyacrylate), covered with a cellulose thin layer resting on an impermeable non woven layer.

Liquids leak prevention.

- 5 layers: absorbing 200 to 300 times its mass in water.
- Neutralize the odours.
- Antimicrobial protection.
- Ideal for dentals, cages, surgery tables and puppy schools.

Sold in packs of 25 or 30 under pads.



EXTREMELY
ABSORBENT !

Code	Pack.	Size	Absorption capacity
327050	x 30	40 x 60 cm	800 mL
327052	x 30	60 x 60 cm	1 000 mL
327054	x 25	60 x 90 cm	1 900 mL
327056	x 25	60 x 90 cm	2 500 mL



DRYSOFT

ADVANCED VETCARE, VIC

Greta, Practice Manager

Genia puppy pads are really good, we never had any leak with, as they are highly absorbent. We like the fact they come in pack of 25.

SUCCESS VET HOSPITAL, WA

Brie, Practice Manager

The Drysoft are very absorbent; good quality.

TINANA VET SURGERY, QLD

Trinity, Senior Head Nurse

Drysoft Ultra puppy pads are the best on the market, excellent quality, no leakage as they absorb a huge amounts of fluids. We won't use anything else. Fantastic product.



TABLET CUTTER & CRUSHER

Our new tablet crusher is portable and easy to clean. Simply put the tablet in the device and screw the top to crush it. Spare tablets can be stored in the middle compartment.

Sold individually.




Code	Size	Diameter
227086	64 mm	40 mm

BOTTLE TEATS

Non-perforated, natural latex teats for kittens, puppies and lambs. Fit on bottle necks. Certified 93/11/CE.

Sold individually.

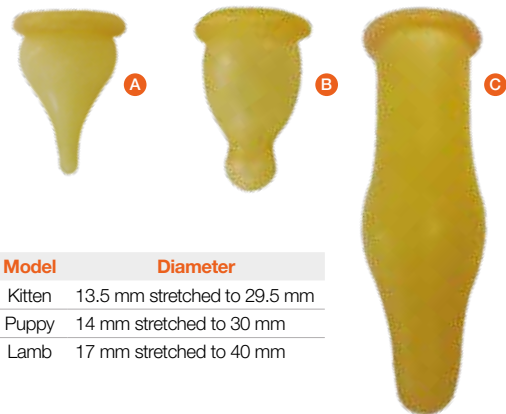

	Code	Model	Diameter
A	227021	Kitten	13.5 mm stretched to 29.5 mm
B	227022	Puppy	14 mm stretched to 30 mm
C	227023	Lamb	17 mm stretched to 40 mm

BLUE PET PILLER

Pill popper with a soft atraumatic rubber tip end. Can be used with water to wash down the tablet.

- Flexy soft tip with good grip.
- Can be used with one hand.
- Add water into the syringe to facilitate the swallowing of the tablet.

Sold individually or in boxes of 12 pet pillers.



Code	Description
227230	Pill popper

3 EASY STEPS:

1. Draw water into the pet piller
2. Add the pill into the rubber tip
3. Flush the tablet down the animal throat



PET PILLER

KINGS MEADOWS VET CLINIC, TAS

Kay, Practice Manager

We like the size of the Pet Piller. We find it easy to use.

MAIN RIDGE VET, VIC

Jacqui, Practice Manager

It's easy to use, I mix it with water.

KYOGLE VET CLINIC, VIC

Rebecca, Practice Manager

It's convenient to use it with water. Really easy to use.

DOG ROCK VET CLINIC, WA

Sam, Practice Manager

Genia Pet Piller is quite good! We like the fact we can use it with water; it encourages the patient to swallow, instead of spitting out the tablet.

MELBOURNE CAT VET, VIC

Samara, Practice Manager

We like the fact Genia Pet Piller are smaller than regular ones, making it easier to use.

HEYWOOD VET CLINIC, VIC



Dr Tracy, Vet

The pet piller is really user-friendly. Pet owners find it easy to use, and don't come back complaining about it. Thanks to the water feature, it helps the automatic swallow reflex, which makes it less stressful for the patient.





THE CHOICE DEPENDS ON THE VETERINARIAN'S CLINICAL ASSESSMENT

TYPES OF WOUND	DRESSINGS	DRESSING CHANGE
 <p>Wound in debridement stage or in granulation stage. Low/medium exudate wound without loss of substance.</p>	<p>HYDRO</p>	<p>Every three days on average. Change if discolored or swollen.</p>
 <p>Wound from debridement to epidermization. Medium/high exudate wound aural hematoma.</p>	<p>SPONGI</p> <p><i>Spongy side on the wound</i></p>	<p>When saturated or every five days.</p>

ANIPLAST® HYDRO

Healing dressing. Second skin.

- Gelatine-free hydrocolloid layer applied to a waterproof polyurethane film.
- Maintains a moist and warm environment at the surface of the wound: promotes granulation and epithelialization along with wound self-debridement.
- Impermeable to bacteria.
- Sticky side on the wound.

Code	Pack.	Dimensions
327014	x 10	10 x 10 cm
327016	x 5	20 x 20 cm



ANIPLAST® SPONGI

Super absorbent dressing. Non-adhesive, semi-occlusive foam dressing on a thick water-resistant polyurethane film.

- To treat chronic and exuding wounds. Maintains a moist and warm wound environment, while limiting maceration and infection.
- Non-adherent: the dressing is suitable to use on fragile skin.
- Facilitate granulation and epithelialization.
- White foam on the wound.

Code	Pack.	Dimensions
327000	x 10	10 x 10 cm
327002	x 3	20 x 20 cm



ABSORBS UP TO 50 TIMES ITS WEIGHT | IDEAL FOR AURAL HEMATOMA |



MONOFIL



Synthetic absorbable monofilament thread (purple).

Composition: polydioxanone (PDO/PDS).

- Ideal for smooth passage through tissue. Soft tissues and hollow organs. Slow healing tissues. GI and urological surgeries. Subcutaneous or intracutaneous suturing.
- Sterilisation method: EO
- Tolerance: no wick effect.
- Minimal sawing effect.
- Excellent knot security.
- Loss of tensile strength: 25% after 30 days.
- Complete absorption by hydrolysis: 180 to 240 days.
- Instructions: should be stored in a cool place out of direct sunlight.
- Use within **4 months after opening**.

Sold individually.

DISPENSER CASE		
USP	15 m	25 m
3/0	340904	340905
2/0	340906	340907
0	-	340909
1	340910	340911

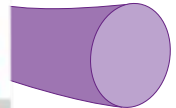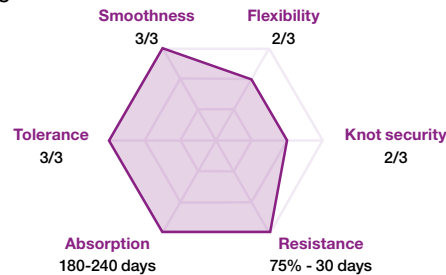

GLYCOFIL®

Synthetic absorbable braided thread (purple).



Composition: polyglycolic acid (PGA) with biocompatible coating.

Coating: polycaprolactone and calcium stearate.

- Ideal for keeping knots secured and easy handling. Muscles, subcutaneous tissues, vascular surgery, ophthalmic surgery, uterus and GI surgeries.
- Sterilisation method: EO
- Smoothness: easy passage through tissues. Minimal sawing effect.
- Loss of tensile strength: 50% after 21 days.
- Complete absorption by hydrolysis: 60 to 90 days.
- Instructions: should be stored in a cool place out of direct sunlight.
- Use within **4 months after opening**.

Sold individually.

DISPENSER CASE	
USP	25 m
3/0	340811
2/0	340813
0	340815
1	340817

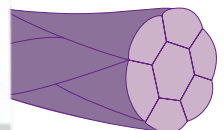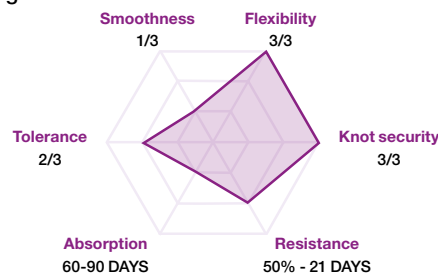

FILAMID

Synthetic non-absorbable pseudomonofilament and coated suture made of polyamide.


Black thread: easier to find the suture when removing it.

Sterilised with Ethylene Oxide (EO).

- Skin suture. Ideal for smooth passage through the skin.
- Sterilisation method: EO

Sold individually.

DISPENSER CASE	
USP	100 m
3/0	340570
2/0 (thinner)	340568
2/0	340566
0	340564

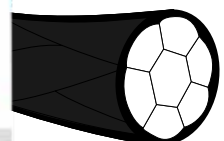



QUALITY CONTROL

STERILIZED IN POUCH

One-step sterilisation process.

REINFORCED

More resistant.

THICK AND WATERPROOF

Minimal hygrometry.




ATTACHED CAP

Avoid losing it.

STORED in an aluminium pouch for minimal hygrometry.

COMPACT

Fit in a pocket.

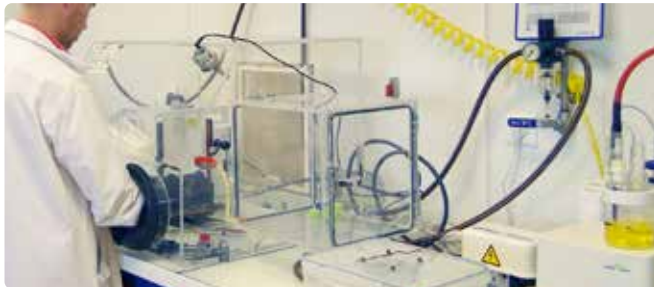

Do not use if external vacuum packaging is pierced.



OK FOR USE



DO NOT USE

TESTIMONIALS

SUTURE

THE RIVERSIDE VET, TAS


Dr Susan, Vet

The suture cassettes are value for money. The thread doesn't have too much memory and stays flat, which is easy to handle.

SUTURE NEEDLES

Stainless steel suture needles, curved, 3/8 circle. Spring double eyes. Triangular or round point.


Sold in boxes of 12 suture needles.

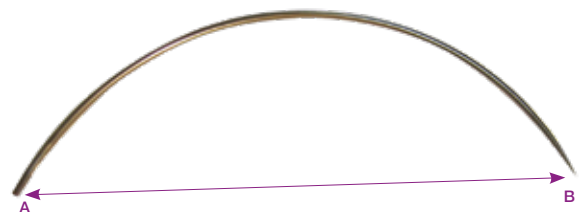

SUTURE-RACK

Suture rack made to hold 3 Génia cassettes. Weighted. Ideal to store all your suture material.

Sold individually.

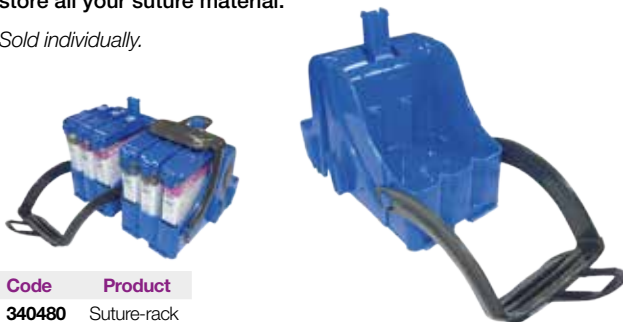

Code **Product**
340480 Suture-rack

Code x 12	Model	Type	Length	Chord (distance A-B)
157723	▲ Triangular point	Spring	45 mm	35 mm
157725	▲ Triangular point	Spring	30 mm	24 mm
157726	▲ Triangular point	Spring	20 mm	16 mm
157742	○ Round point	Spring	45 mm	35 mm
157744	○ Round point	Spring	30 mm	24 mm
157745	○ Round point	Spring	20 mm	16 mm

▲ Standard cutting needle: to be used for skin, muscle fascia.
○ Round bodied taper needle: to be used for hollow organs, blood vessels.




DRAPVET

Single-use surgical drapes, non-woven, sterile.

- Convenient and economical.
- One smooth waterproof layer to prevent strike-through.
- One non-woven, comfortable and absorbent layer. Soft material.
- Made of polypropylene and polyethylene.
- 53g/m².

Sold in boxes of 6, 12 or 24 drapes. Singly sterilised and wrapped.

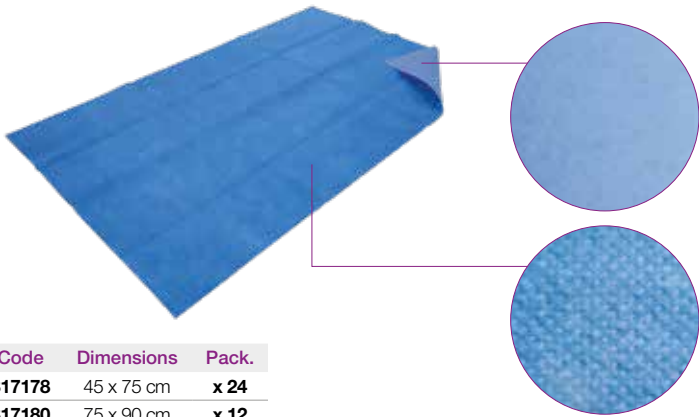

Code	Dimensions	Pack.
317178	45 x 75 cm	x 24
317180	75 x 90 cm	x 12
317182	125 x 125 cm	x 12
317184	200 x 300 cm	x 6



STAPVET

Sterile disposable 35W skin stapler.

- 35 stainless steel staples (7 mm x 4 mm).
- Angled head for good visibility.
- Ergonomic handle.
- Easy and precise staple placement.
- Non bouncing.

Sold individually. Sterile.




Code	Model
2278601	35W stapler

STAPLE REMOVER

Disposable skin staple remover.

- The round and curved ends are designed to preserve the skin and to easily fit under staples.
- Ergonomic.

Sold individually. Sterile.

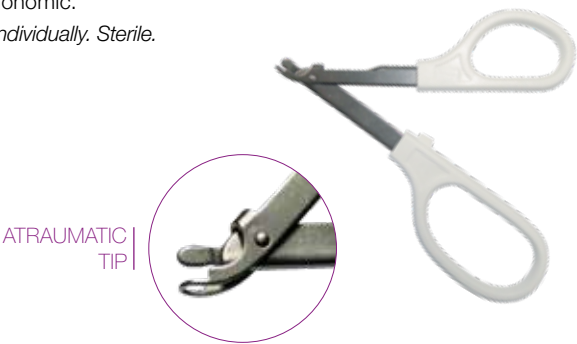

Code	Model
2278621	Staple remover



HANDY CANIS


Adjustable and washable harness / sling for hind legs. Soft and comfortable for your patients thanks to its stretchable heavy-duty neoprene material. Very easy to adjust with a wide Velcro closure secured by a second fastening system.

- Ideal to help animals going up and down stairs or getting in and out of your car.
- Recommended for hydrotherapy.
- Help patients suffering from arthritis, degenerative myelopathy or paralysis, walking.
- Can be used to lift heavy dogs, instead of using towels.
- Adjustable lifting straps for the owner comfort, as well as a padded handle for a snug grip. **A**
- Entirely washable.

All straps can be adjusted to accommodate breeds of varying heights and weights. Sold individually.



Code	Size	D Abdomen circumference min. and max.	E Thigh circumference min. and max.	Breed
227250	XS	30 - 33 cm	7,5 - 17 cm	Yorkshire, Chihuahua
227252	S	33 - 36 cm	15 - 23 cm	Dachshund
227254	M	36 - 50 cm	22 - 32 cm	Basset Hound
227256	L	50 - 63 cm	31 - 36 cm	Golden Retriever, German Shepherd
227258	XL	63 - 81 cm	37 - 53 cm	St. Bernard, Newfoundland

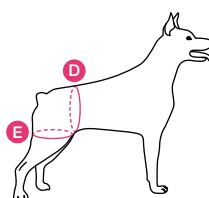

D) Abdomen circumference at groin.
E) Thigh circumference at groin.



NEOPRENE
MACHINE WASHABLE

CANIS-MOBILE®

Kart used in case of mobility loss to support dogs' hind leg(s). Easy to set up and completely adjustable.

- Welded, galvanised metal and aluminium chassis.
- Adjustable height, length and straps.
- Soft rubber wheels.
- Removable neoprene harness to ease washing.

Indications:

- Paralysis or paresis of dogs' hind legs: great for German Shepherds, Dachshunds, Labradors.
- Functional rehabilitation after diagnosed herniated disc.
- Temporary immobilisation post cruciate ligament surgery or injury, cast, splint.


Supplied with a hind legs harness (Handy Canis).

CANIS-MOBILE® measurements

- A) Thorax circumference behind the front legs.
B) Length between the front and back legs.
C) Tail height.

HANDY CANIS measurements

- D) Abdomen circumference at the groin.
E) Thigh circumference at the groin.

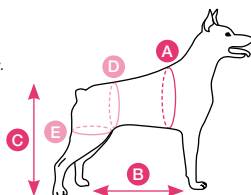

Help with Canis-Mobile® selection: measure A, B and C and determine the kart size (small, medium or large). Then select the compatible hind legs harness size (D and E).

Code SET	HANDY CANIS size	HANDY CANIS measurements		CANIS MOBILE® size	CANIS-MOBILE® measurements			Max. Weight
		D Abdomen circumference min. and max.	E Thigh circumference min. and max.		A Thorax circumference	B Length min. and max.	C Tail height min. and max.	
260571	S	33 - 36 cm	15 - 23 cm	Small	41 cm maxi	16 - 26 cm	17 - 27 cm	10 kg
260575	M	36 - 50 cm	22 - 32 cm	Medium	66 cm maxi	25 - 35 cm	24 - 38 cm	30 kg
260580	L	50 - 63 cm	31 - 36 cm	Large	90 cm maxi	35 - 50 cm	41 - 56 cm	50 kg
260581	XL	63 - 81 cm	37 - 53 cm					





EXAMED

Latex pre-powdered disposable examination gloves. Lightweight for improved dexterity and reduced hand fatigue.

- Ambidextrous, strong and comfortable gloves.
- Extremely stretchable (+650%): no more tearing !
- Rolled edges.
- Non sterile.

Sold in boxes of 100 gloves.



LIGHTLY POWDERED
EXAMINATION GLOVES

LATEX |



Code x 100	Size
307760	6/7 S
307761	7/8 M
307762	8/9 L
307763	9/10 XL

TESTIMONIALS

GLOVES

ARMSTRONG CREEK VET PRACTICE, VIC

Rachel, Vet Nurse

We never had any problem with Genia gloves. They are comfortable and smooth, without leaving your hands dry.

NEW

EXAMED POWDER-FREE

Powder-free disposable examination gloves made of latex.

- Ambidextrous, strong and comfortable gloves.
- Good stretchability > 650% elongation.
- Double dipped to reduce porosity.
- Rolled edges.
- Non-sterile.

Sold in boxes of 100 gloves.



POWDER-FREE |

EXAMINATION
GLOVES

LATEX |



Code x 100	Size
307745 ●	6/7
307746 ●	7/8
307747 ●	8/9
307748 ●	9/10

NEW

VETRILE

Single-use synthetic examination gloves, 100% nitrile, powder-free, non-sterile.

- Non spermicidal effect.
- Ambidextrous.
- Beaded cuff increasing the gloves strength. Easier to put on.
- Limited allergic risk: no latex protein (type IV allergies due to chemical accelerators possible).
- Good stretchability (400%) and resistance to perforations, as well as chemicals.
- Antimicrobial barrier.

Sold in boxes of 100 gloves.



NITRILE
100%

WITHOUT
LATEX |

LIMITED
ALLERGENICITY |



Code x 100	Size
307769 ●	6/7
307770 ●	7/8
307771 ●	8/9
307772 ●	9/10





NEW

SURGITEK

Sterile latex surgical gloves, lightly powdered.

Ideal thickness for strong but sensitive gloves with excellent stretchability to add extra comfort and durability. Smooth and microtextured on the palm to optimise the grip on instruments.

- High quality natural latex with low protein content to limit allergenicity.
- Ultra sensitive for an optimal touch.
- Micro-textured surface for a better grip.
- Anatomically-designed fingers and thumb on the inside.
- Rolled edges and lightly powdered for easier removal.

Sold in boxes of 50 pairs. Singly sterilised and wrapped.



LIGHTLY POWDERED |

SURGICAL
GLOVES |

LATEX |



Code x 50	Size
307698	6
307700	6.5
307702	7
307704	7.5
307706	8
307708	8.5

SURGITEK POWDER-FREE

Anatomically-designed powder-free sterile surgical gloves with long cuffs.

- High quality natural latex with low protein content to limit allergenicity.
- Powder-free.
- Ultra sensitive for an optimal touch.
- Micro-textured surface for a better grip.
- High breaking strength: Median breaking strength 9 N minimum.
- Rolled edges for easier removal.

Sold in boxes of 50 pairs of gloves. Singly sterilised and wrapped.



POWDER-FREE |

SURGICAL
GLOVES |

LATEX |



Code x 100	Size
307710 ●	6
307712 ●	6.5
307714 ●	7
307716 ●	7.5
307718 ●	8
307720 ●	8.5







Our range of admin sets are known for the quality of their tubing, as well as their multiple injection ports.
All compatible with most IV pumps on the market.




LUER-LOCK HUB

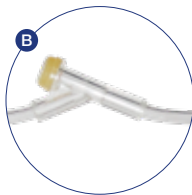

INJECTION PORT

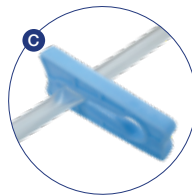

CLAMP

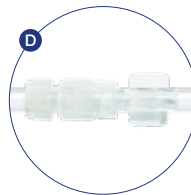

LUER-LOCK CONNECTION

PERFUVET PLUG

RIVERTON ROSSMOYNE VET HOSPITAL, WA
Stacey, Practice Manager

The vets are happy with the size & length of the Perfuvet Plug and like the extra ports.

PERFUVET DUO

PANANIA VET CLINIC, NSW
Leanne, Practice Manager

Better quality than our previous product. Not pulling apart at connection for portals. 2 portals as compared to some other giving set.

PERFUVET MICRO

GARDEN CITIES VET CLINIC, WA
Dr Simon, Vet

I like the Perfuvet Micro giving set because we don't have any IV pump. It is lovely to have some degree of metered control of the flow rate. We do make sure to hang the I/V bag as high as possible to maximise the accuracy. Works well in small to medium dogs and cats.

BYFORD VET HOSPITAL, WA
Brooke, Practice Manager

The Perfuvet Micro is convenient; we use it for our anesthetic and hospitalised patients.

TESTIMONIALS

PERFUVET® DUO

Premium IV admin set for small animals with two injection ports. Compatible with most IV pumps. Premium quality tubing to prevent kinking.

- Calibration: **20 drops/ml.**
- 200 cm Admin set: PVC tubing.
- Luer-Lock.
- Two injection site ports at 10 cm and 110 cm from the catheter.
- Convenient clear and soft drip chamber.
- On/off clamp.
- Sharp piercing spike.
- Flow control roller clamp.

Sold in boxes of 12 IV sets.

2 INJECTION PORTS

Code	Product
157395	PERFUVET® DUO

PERFUVET® PLUG


IV admin set for small animals with two injection ports. PERFUVET® PLUG is a 2 in 1 admin set with an extension set. Disconnect the extension line at 25 cm from the catheter. Compatible with most IV pumps. Premium quality tubing to prevent kinking.

- Calibration: **20 drops/ml.**
- 200 cm Admin set: PVC tubing.
- Luer-Lock.
- Two injection site ports at 10 cm and 110 cm from the catheter.
- Convenient clear and soft drip chamber.
- On/off clamp.
- Sharp piercing spike.
- Flow control roller clamp.

Sold in boxes of 12 IV sets.

25 CM EXTENSION SET INCLUDED: 2 IN 1 !

Code	Product
157397	PERFUVET® PLUG

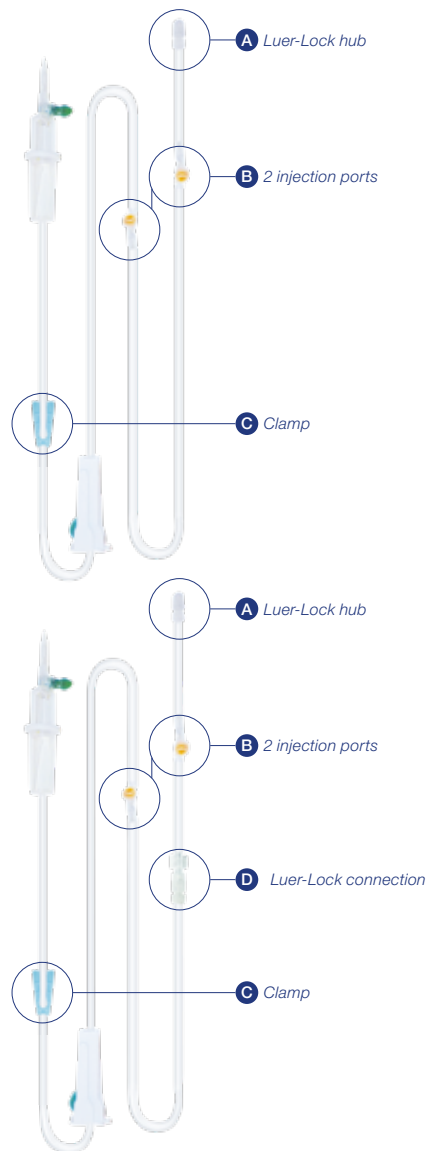



PERFUVET® MICRO


IV admin set for small animals with two injection ports. The dial-type flow controller allows an easy and accurate setting of drip rates: ideal for paediatrics, or when you run out of IV pumps. Manually set up the ml/h.

PERFUVET® MICRO is also a 2 in 1 as it includes a 25 cm extension set.

- Calibration: **60 drops/ml**.
- 200 cm Admin set: PVC tubing.
- Luer-Lock.
- Convenient clear and soft drip chamber.
- On/off clamp.
- Sharp piercing spike.
- Great to use for CRI's (Constant Rate Infusions).

Sold in boxes of 12 IV sets.

Code	Product
157398	PERFUVET® MICRO

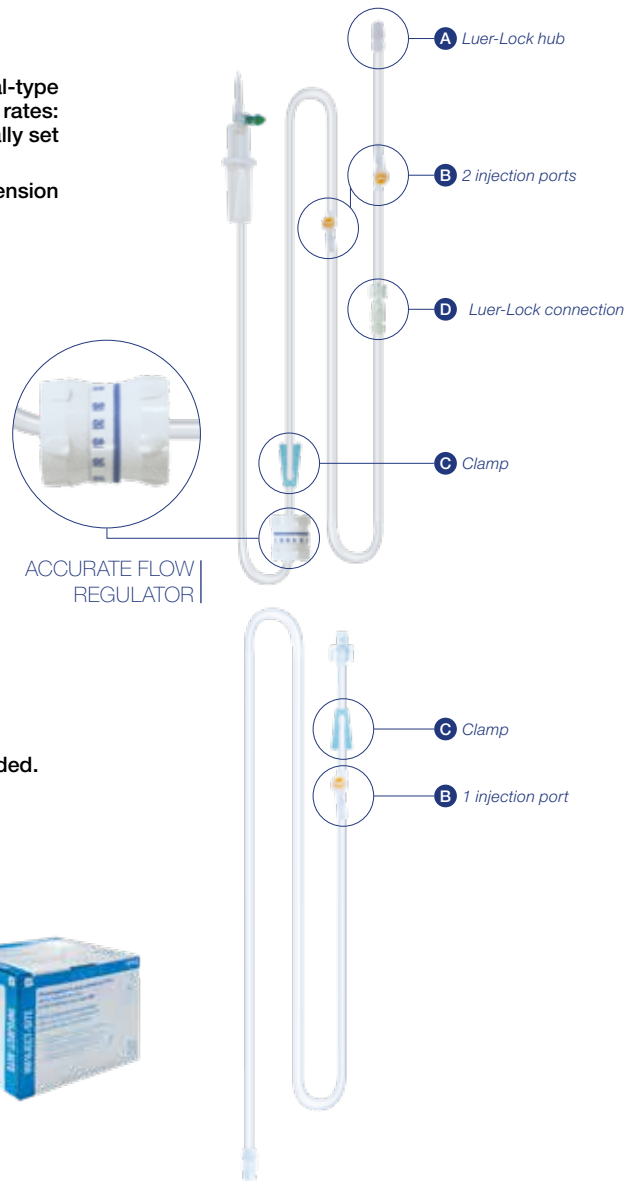

INFUJECT-SITE

Extension sets (26,5 cm and 100 cm) with injection port included.

- PVC extensions with Luer-Lock hub.
- Injection port located at 5 cm from the catheter.
- On/off clamp.
- Available in 2 lengths.

Sold in boxes of 25 sets.

Code	Diameter		Length	Tube capacity
	Int.	Ext.		
157399	2.9 mm	4.1 mm	26.5 cm	2.08 mL
157400	2.9 mm	4.1 mm	100 cm	6.72 mL



K-STAR

Star-shaped device, to use on admin sets to avoid kinking. It helps secure the catheter, and facilitates fluid flow. Quick and easy, it fits on standard tubings. Can also be used as a J Loop connector.

Sold in boxes of 50 K-stars.



Code	Pack.	Product
191052	x 50	K-STAR - Anti kinking device

INJECTION PORT

Small see-through injection ports with luer-lock and wide injection platform.

Sold in boxes of 12 injection ports. Singly sterilised and wrapped.



Code	Description	Pack.
157330	Male cap	x 12

TESTIMONIALS

K-STAR

CAPE CREATURES VETERINARY HOSPITAL, WA

Genene, Head Vet Nurse

The K-star is fantastic to keep IV drops and avoid tubing kinking; really easy to use AND re-usable.

ABERDEEN VET CLINIC, VIC

Sharron, Practice Manager

We love the K-star! It's ideal to stop the IV line from kinking.

GARDEN CITIES VET CLINIC, WA

Dr Simon, Vet

The K-star is fantastic, really simple. Nice tool.



PERFUVET® TWIST

Let's take away some of your day-to-day worries with our new PERFUVET® TWIST. This new device prevents patients from getting tangled up in the tubing. Composed of a reusable rotating base placed on the cage ceiling and a spiral extension set, this is the best combination for patients to move freely and prevent fluid flow obstruction. Why do you love it?

- When patients walk in circles in the cage the spiral tubing will follow their movements, thanks to the rotating base.
- Tubing will not tangle around the animal and will be less inclined to kink thanks to its spiral shape.
- It is ideal for hospitalized patients, overnight IV therapy, and patients with neurological disorders.
- It will improve safety and save time.



Ready to use kit includes:
- One grey base with an adhesive backing for permanent attachment to the top of the cage.
- Two blue shells connecting the extension set to the base.
- Three sterile 3 meters spiral extension sets with collar ring.

Sold as a kit including one base, two blue shells and three sterile spiral extension sets.



Code	Description
157450	Perfuvet® Twist kit (1 base, 2 blue shells and 3 spiral tubings)
157451	3 spiral tubings
157452	Set of 2 blue shells

IDEAL FOR YOUR:

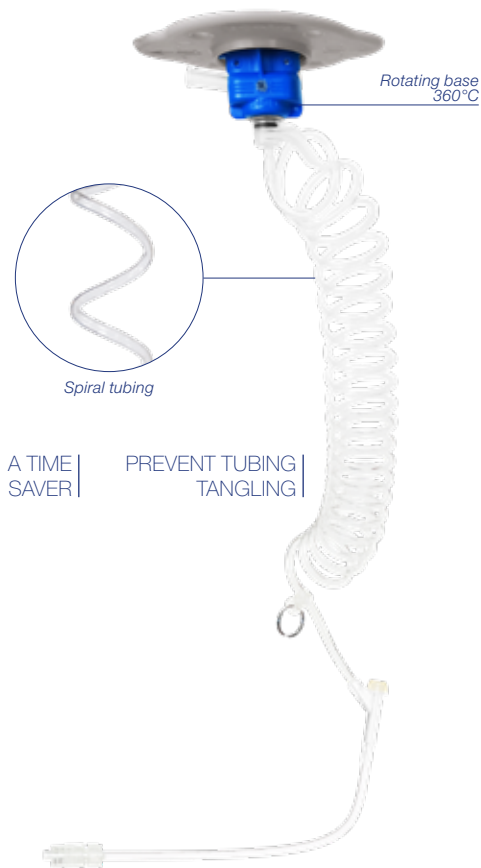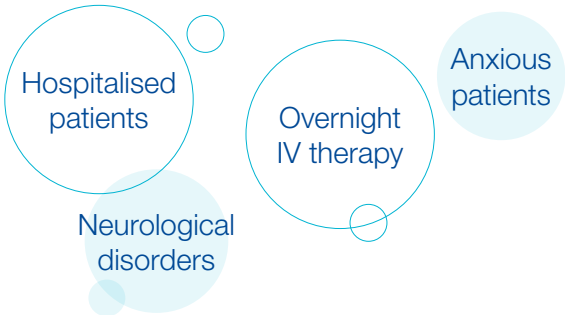

A TIME SAVER | PREVENT TUBING TANGLING

PERFUVET TWIST

WYNYARD VET CLINIC, TAS

Jane, Vet Nurse

The Perfuvet Twist is brilliant! We used it with a hospitalised dog who didn't get tangled in the cage.

COMPANION ANIMAL PRACTICE, NSW

Dylan, Head Vet Nurse

We have definitely been using the Perfuvet twist. They are a great product, each time we use them they are easy to install and the animals are about to move around the cage freely.


ABERDEEN VET CLINIC, VIC

Sharron, Practice Manager

We use the Perfuvet Twist for our hospitalised patients. We find it much easier than regular admin sets, as it's not getting tangled around the patient's legs.



1

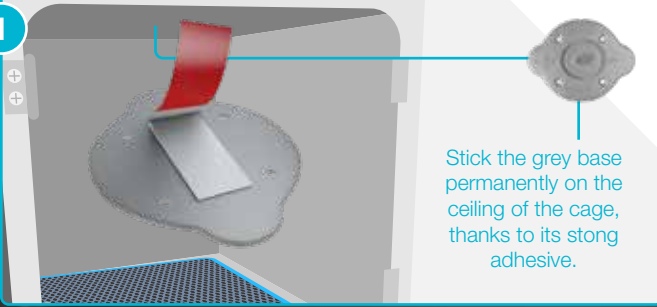


2




3




4

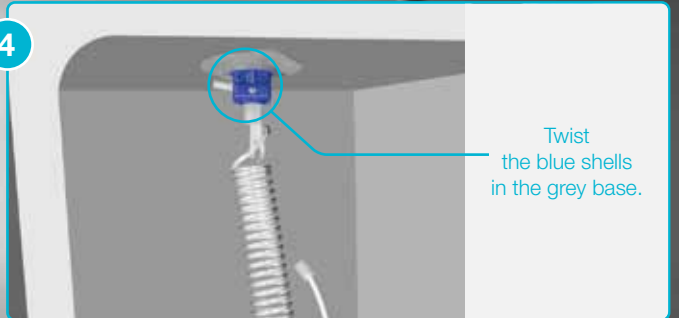

5

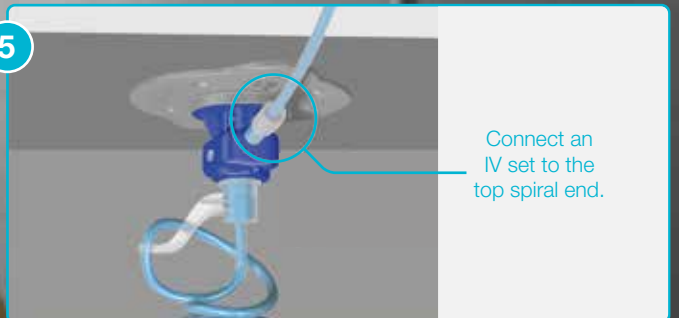

6

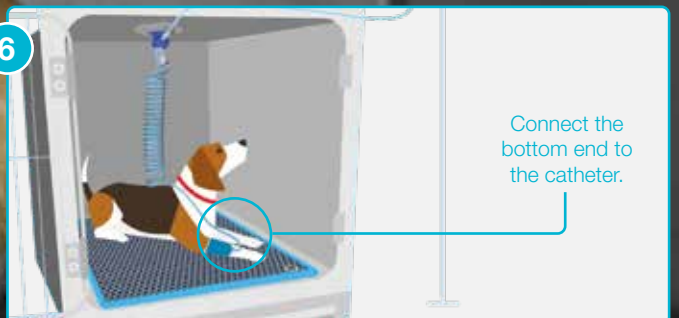





ANIVEN®

FEP intravenous catheter featuring flexible wings for easy placement. Ultra-sharp tri-beveled needle.

- Long triple-bevel designed to minimise patient discomfort and tissue trauma **A**.
- A silicone layer coating the needle and catheter results in a smooth transition from needle to catheter during insertion.
- Silicone-coated Teflon Fluoroethylene Propylene (FEP) catheter to prevent kinking during catheter insertion **B**.
- Transparent viewing chamber for immediate visualisation of blood reflux.
- Porous membrane allows a quick blood rise into the viewing chamber **C**.
- Ergonomic 'finger rest' for a smooth and controlled placement **D**.
- Movable and flexible wings with two grooves to comfortably fit animal's limb **E**.
- Catheter can be stitched in place, thanks to its wing holes.
- Luer-Lock with the catheter colour code **F**.

Sold in boxes of 12 catheters. Singly sterilised and wrapped.

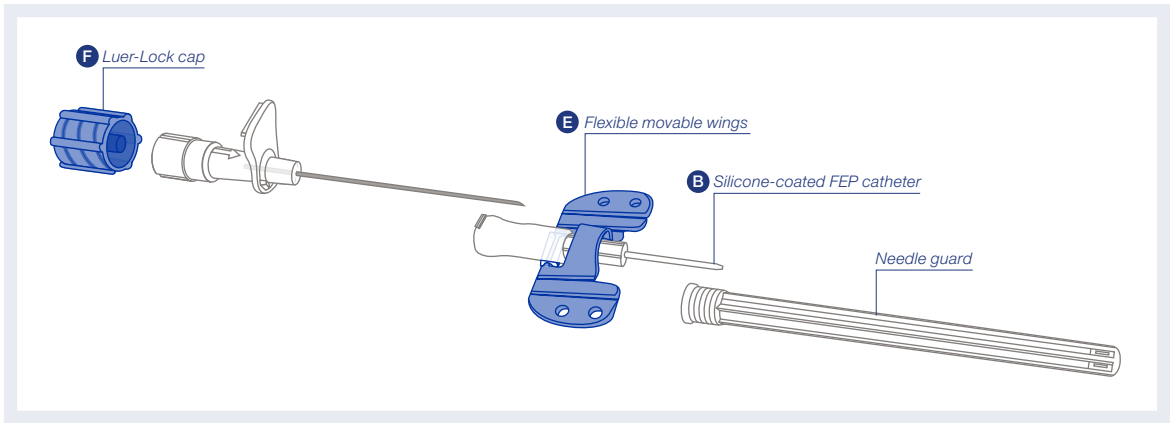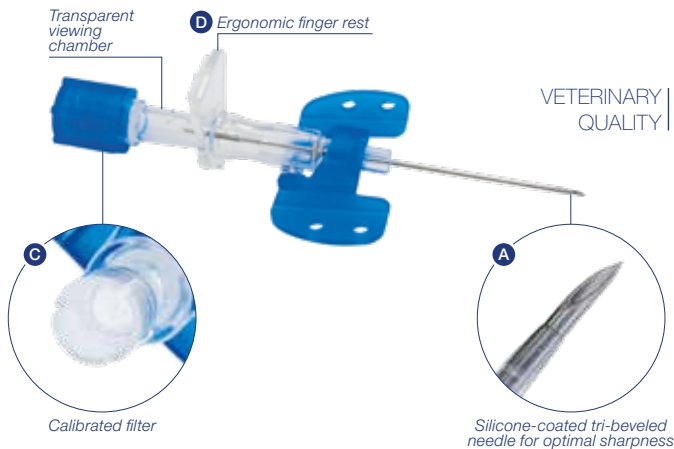


TESTIMONIALS

ANIVEN

ADVANCED VET NURSE TRAINING, WA

Samantha, Owner

I tried out a new style of IV catheter called Aniven from Genia Australia. I think these will be great for those patients who roll around in their cage and pull their catheter out, the wings will be much easier to tape in than other types of catheters.

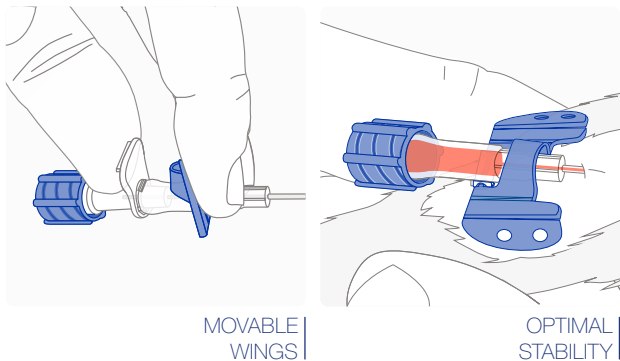

Code x 12	Length	Size (G)	Ext. diameter
157850	19 mm	24G x 3/4"	0.7 mm
157852	25 mm	22G x 1"	0.9 mm
157854	33 mm	20G x 1 1/4"	1.1 mm
157856	45 mm	18G x 1 3/4"	1.3 mm








ANIVASC®

FEP intravenous catheter with fixed wings and a tri-beveled needle.


- Long tri-beveled needle for great pricking, efficient and painless insertion **A**.
- Optimal glide thanks to the silicone covering the needle and the catheter.
- Silicone-coated Fluorinated Ethylene Propylene (FEP) catheter to prevent from kinking during catheter insertion **B**.
- Transparent viewing chamber for immediate visualisation of blood reflux.
- Porous membrane allowing a quick blood rise into the viewing chamber **C**.
- Ergonomic 'finger rest' for a smooth and controlled placement **D**.
- Fixed colour-coded polypropylene wings **E**.
- Catheter can be stitched on, thanks to the holes on the wings.
- Luer-Lock cap **F**.

Sold in boxes of 12 catheters. Singly sterilised and wrapped.

anivasc
veterinary excellence

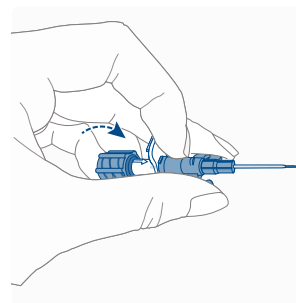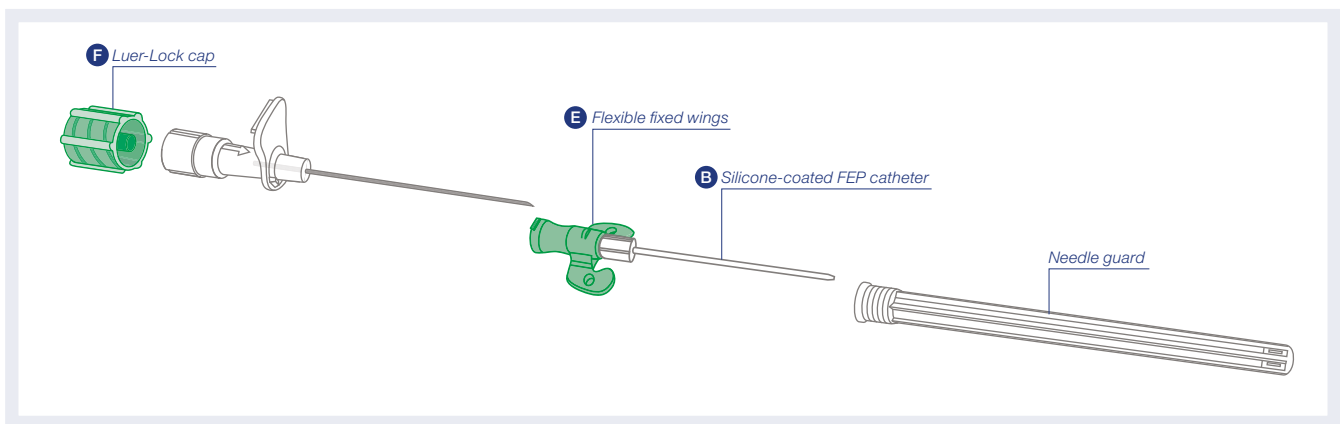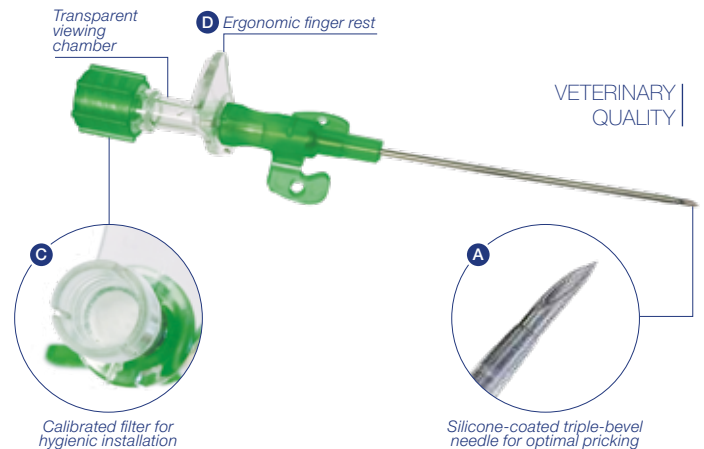

COMFORTABLE
GRIP

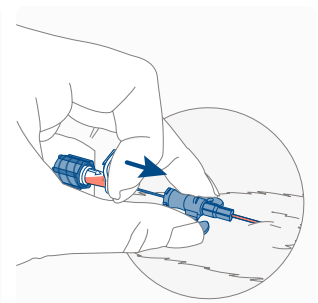

OPTIMAL
INSERTION

Code x 12	Dimensions	Size (G)
157900	0,7 x 19 mm	24G
157902	0,9 x 25 mm	22G
157904	1,1 x 33 mm	20G
157906	1,3 x 45 mm	18G
157864	1,6 x 105 mm	16G
157868	2,1 x 80 mm	14G
157870	2,1 x 105 mm	14G





DISPOJECT® SAFETY INSULIN SYRINGE

40 IU insulin syringe of 1 mL with a 29G x 1/2" (13 x 0.3 mm) ultra-fine needle.

- Sterile, disposable, latex-free.
- Innovative safe system to secure the needle after use in three easy steps:
 - After the injection, pull the plastic guard up
 - Twist it to lock it
 - Snap the top off
- Save space in your sharp container !

Sold in boxes of 100 syringes. Singly sterilised and wrapped.


1. After the injection, pull the plastic guard up.



2. Twist it to lock it.



3. Snap the top off.

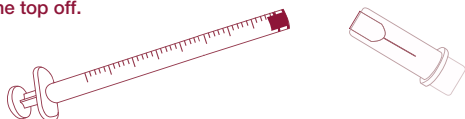

Code x 100	Size	Graduation
107136	1 mL	40 IU



**TWIST IT
TO LOCK IT !**



PLEASE CONTACT US
TO GET YOUR
FREE INFORMATIVE POSTER
FOR YOUR CLIENTS






PRIMAJECT® NEEDLES

Stainless steel disposable hypodermic needles.

- Tri-beveled and silicone-coated cannula.
- 6% taper, Luer-Lock.
- Epoxy resin hub.
- Sterilised with ethylene oxide.
- Latex-free, DEHP-free.

Sold in boxes of 100 needles.

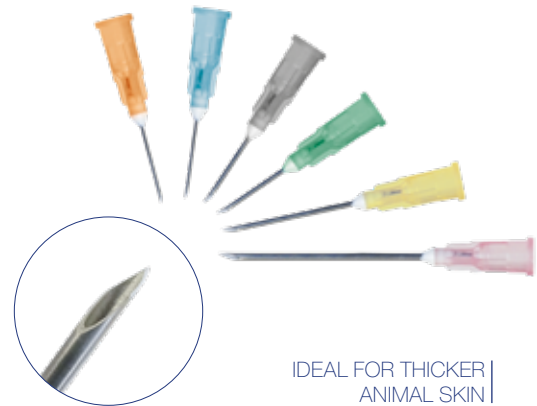

IDEAL FOR THICKER |
ANIMAL SKIN |

MADE FOR |
VETS |

Code x 100	Size	Size (G)
147602	0,5 x 16 mm	25G x 5/8"
147604	0,6 x 16 mm	23G x 5/8"
147608	0,6 x 25 mm	23G x 1"
147610	0,7 x 16 mm	22G x 5/8"
147616	0,8 x 16 mm	21G x 5/8"
147618	0,8 x 25 mm	21G x 1"
147624	0,9 x 25 mm	20G x 1"
147634	1,2 x 38 mm	18G x 1 1/2"







UROLENE® CAT

High quality sterile urinary cat catheter featuring a rotating luer-lock end, providing better comfort for the patient.

- **360° rotating collar** - with four stitching holes. Offers more comfort on patient's prepuce, as it follows the cat's movement.
- **For a perfect fit, choose between:**
 - **Closed rounded end with two side holes:** atraumatic, and easier to insert into the urethra.
 - **Open end:** ideal to unblock a cat with urethral obstruction.
- **Luer-Lock hub:** easier for urine collection.
- Made of medical grade elastomer of Polyether Block Amide providing a great balance between rigidity and flexibility. Non-reactive material, ideal to leave in, up to three days.
- **DEHP-free & radiotransparent.**
- Available in **3 sizes: open or closed end and colour coded (blue: 2.5FG, pink: 3FG, green: 4 FG).**

Sold individually or in boxes of 12 catheters. Singly sterilised and wrapped.

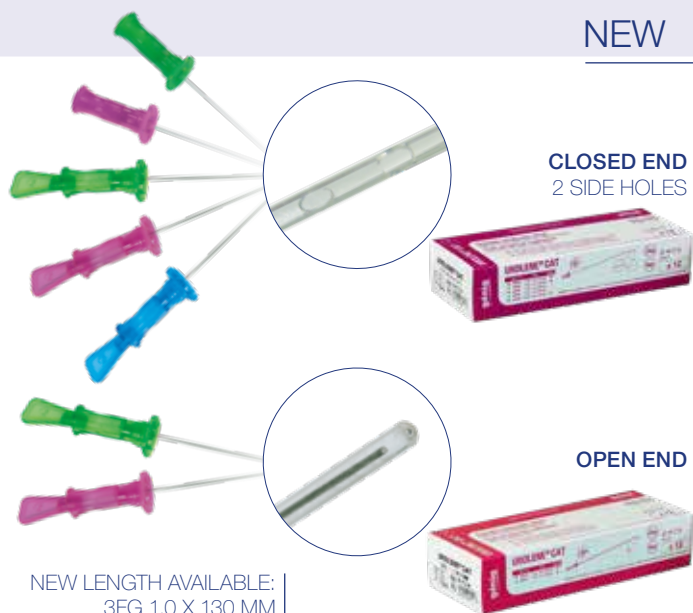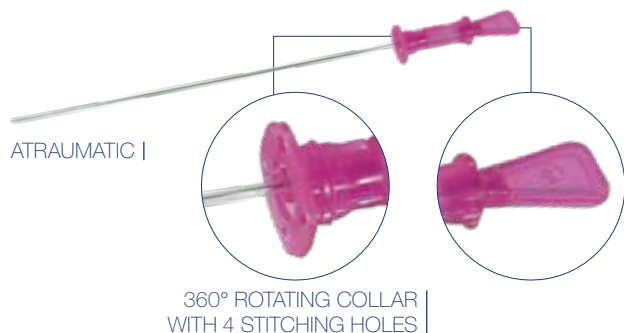



Code Closed end	Code Open end	Model	Size
306898	-	With stylet	2.5 FG 0.8 x 110 mm
306900	306940	With stylet	3 FG 1.0 x 110 mm
306906	306941	With stylet	3 FG 1.0 x 130 mm
306902	306942	With stylet	4 FG 1.3 x 140 mm
306904	-	Without stylet	3 FG 1.0 x 110 mm
306907	-	Without stylet	3 FG 1.0 x 130 mm

UROLENE® DOG

High quality sterile urinary dog catheter.

- **Rounded closed end with two side holes:** atraumatic and easy to insert into the urethra.
- **Excellent balance between flexibility & rigidity.**
- **Luer-Lock hub:** easier for urine collection.
- Made of medical grade Polyamide.
- **DEHP-free & radiotransparent.**
- Available in **5 sizes, colour coded according to the size.**

Sold individually or in boxes of 12 catheters. Singly sterilised and wrapped.

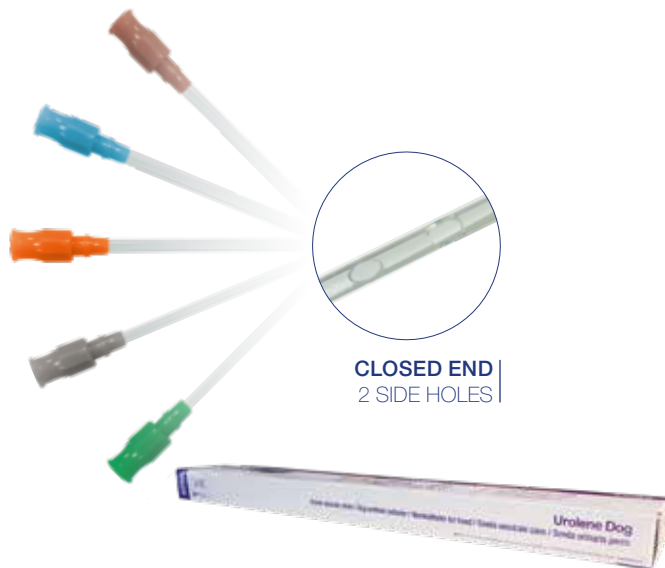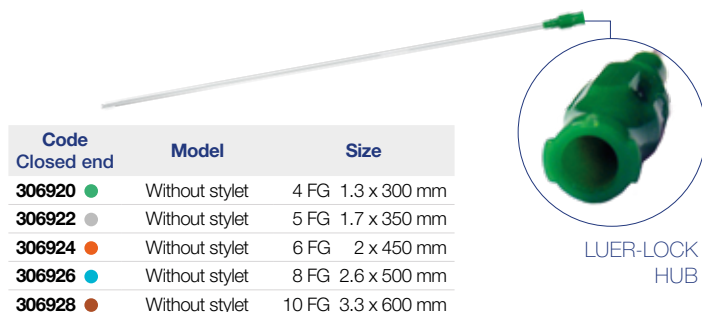

Code Closed end	Model	Size
306920	Without stylet	4 FG 1.3 x 300 mm
306922	Without stylet	5 FG 1.7 x 350 mm
306924	Without stylet	6 FG 2 x 450 mm
306926	Without stylet	8 FG 2.6 x 500 mm
306928	Without stylet	10 FG 3.3 x 600 mm

URINARY CATHETER

ELTHAM VETERINARY PRACTICE, VIC

Amanda, Head Vet Nurse

Urolene catheters are so nice and smooth, compared to others. The material makes them easy to insert into patients.

LISMORE VET CLINIC, NSW

Dr Richard, Vet

We've been using Genia urinary catheters for years. They don't spin too much around at the junction and come in several sizes.

CROYDON VET CLINIC, VIC

Tara, Head Vet Nurse

Genia urinary catheters are easy to work with; the vets like them.

GREENCROSS VETS ARMADALE, WA

Venessa, Practice Manager

Genia urinary catheters are so easy to use! They go through really well. They have as well a wide size reference, to fit them all! Cat catheters: We like the fact we can easily stitch the base in. Dog catheters: They are durable. They're nice and sturdy, without being too stiff.

THE VET GROUP WOLLASTON, VIC

Libbie, Practice Manager

We like the fact we can stitch the urinary catheters. They are really comfortable for the patients.



SURGICAL INSTRUMENTS



METZENBAUM SCISSORS

Surgical quality stainless steel Metzenbaum scissors with shine finish.



	Code	Description
A	601410	14 cm - Cryocut - Screw - Straight
B	601413	14 cm - Cryocut - Screw - Curved
	601412	18 cm - Cryocut - Screw - Straight
	601414	18 cm - Cryocut - Screw - Curved

MAYO STILLE SCISSORS


Surgical quality stainless steel Mayo stille scissors with shine finish.



	Code	Description
	601420	15 cm - Cryocut - Screw - Straight
	601422	17 cm - Cryocut - Screw - Straight
	601423	15 cm - Cryocut - Screw - Curved
	601424	17 cm - Cryocut - Screw - Curved

MAYO SERRATED SCISSORS


Surgical quality stainless steel Mayo serrated scissors with shine finish.

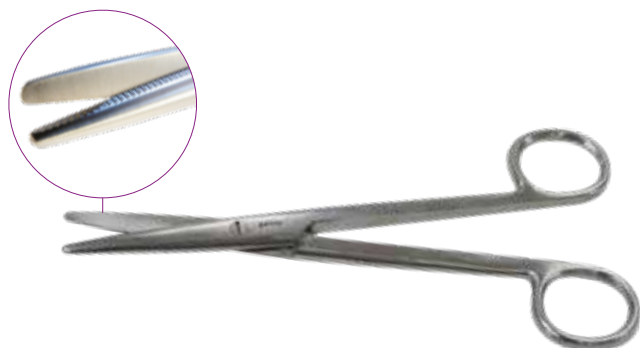

	Code	Description
	601456	17 cm - Surgical - Screw - Straight

SPENCER SCISSORS

Surgical quality stainless steel Spencer scissors with matte finish.

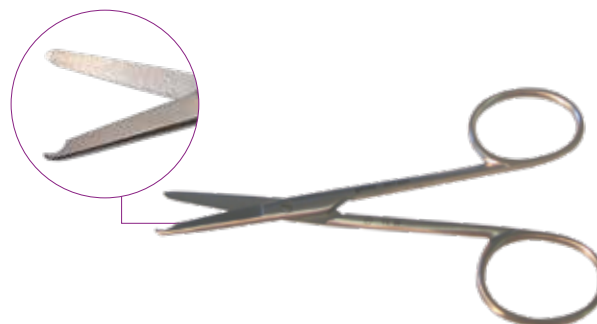

	Code	Description
	601256	11 cm - Surgical - Screw - Straight

IRIDECTOMY SCISSORS

Stainless steel screw-type scissors with satin matte finish.



	Code	Description
	600928	10,5 cm - Surgical - Straight
	600929	10,5 cm - Surgical - Curved



OLSEN HEGAR NEEDLE HOLDER


Stainless steel needle holder with tungsten plates. Matte finish.



Code	Description
607028	14 cm - Tungsten
607030	18 cm - Tungsten

MAYO HEGAR NEEDLE HOLDER


Stainless steel needle holder with tungsten plates. Matte finish.

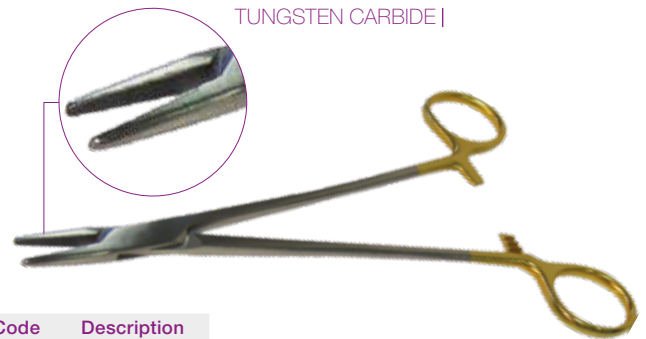

Code	Description
607003	14 cm - Tungsten
607006	18 cm - Tungsten

MATHIEU NEEDLE HOLDER

Stainless steel needle holder with tungsten plates. Matte finish.

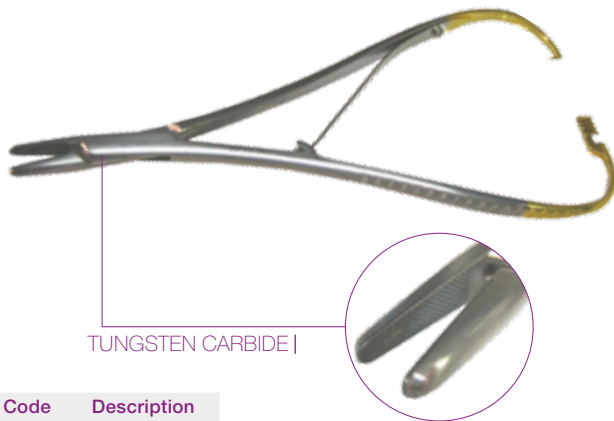

Code	Description
607016	17 cm - Tungsten

GILLIES NEEDLE HOLDER

Stainless steel needle holder. Matte finish.



Code	Description
607022	16 cm - Standard

TOWEL CLAMP

Towel clamps.



Code	Length
650882	9 cm

KILLIAN SPECULUM

Stainless steel 9 cm vaginal speculum for small animals. Set screw adjustment for desired opening.




Code	Length
307980	9 cm



ADSON FORCEPS


Surgical quality stainless steel dissecting forceps with matte finish.

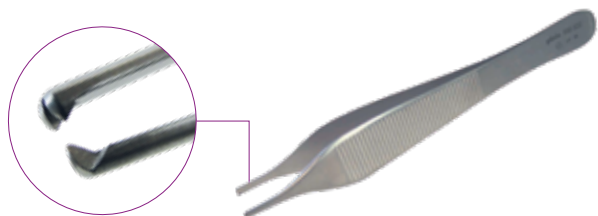

Code	Description
650320	12 cm - Without teeth
650322	12 cm - With teeth

ADSON BROWN FORCEPS


Surgical quality stainless steel dissecting forceps with matte finish.

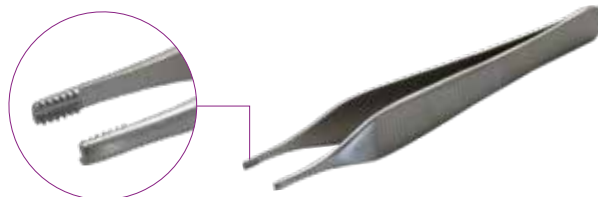

Code	Description
650335	12 cm - With teeth

DEBAKEY FORCEPS

Surgical quality stainless steel dissecting forceps with matte finish.

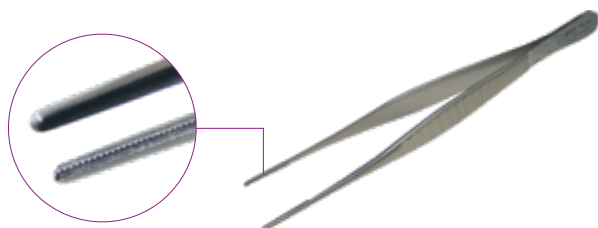

Code	Description
650580	15 cm - Straight - Teeth 1,5 mm

CRILE FORCEPS

Surgical quality stainless steel hemostatic forceps.



Code	Description
657014	14 cm - Without teeth - Straight
657020	14 cm - Without teeth - Curved



HALSTEAD MOSQUITO FORCEPS


Surgical quality stainless steel hemostatic forceps.



Code	Description
657004	12 cm - Without teeth - Straight
657006	12 cm - Without teeth - Curved

ROCHESTER PEAN FORCEPS

Surgical quality stainless steel hemostatic forceps.

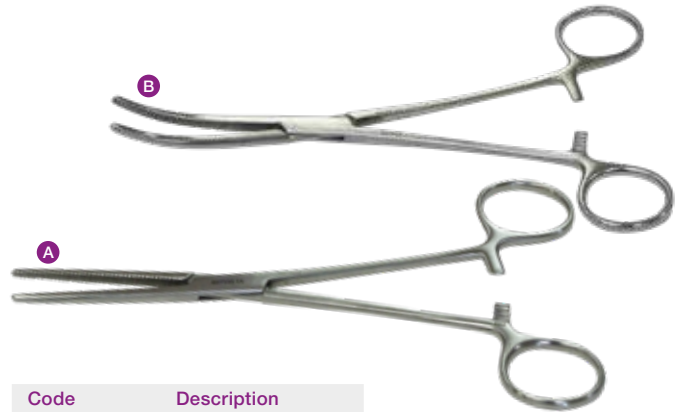

Code	Description
A 657030	18 cm - Without teeth - Straight
B 657032	18 cm - Without teeth - Curved

ALLIS FORCEPS

Surgical quality stainless steel hemostatic forceps.



Code	Description
657000	15,5 cm - With teeth - Straight

HARTMANN FORCEPS - ALLIGATOR

Stainless steel Hartmann foreign body forceps.




Code	Length
307976	8 cm
307975	14 cm

OVARIECTOMY HOOK

Surgical quality ovariectomy hook.

- Round knurled handle for excellent grip.
- Medium ball size (ø 2.8 mm) for optimal grip.

Sold individually.

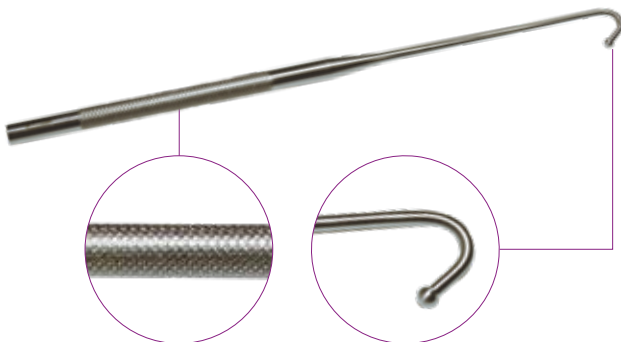

Code	Outside diameter	Length
220360	6 mm	17,5 cm







INTRODUCING MANO MEDICAL: YOUR INNOVATIVE PARTNER IN VETERINARY EQUIPMENT

Mano Medical, now part of the Genia Group, brings over 20 years of expertise in providing high-quality veterinary equipment designed to meet the needs of veterinary practices. Specialising in digital imaging solutions and innovative products, Mano Medical is committed to enhancing patient care and supporting veterinary professionals.

Through Genia Australia, we are pleased to introduce a selection of Mano Medical's cutting-edge products:

Warming System:

Ensure the optimal body temperature of your patients with our advanced carbon technology warming system, specifically designed for veterinary use.

Blood Pressure:

Discover the precision and ease of the Vet BP Doppler, offering reliable and accurate blood pressure monitoring for veterinary practices.

Genia Seal & Cut:

The Genia Seal & Cut is a versatile electrosurgery device that meets a wide range of surgical needs. It combines monopolar and bipolar modes, multiple cutting and coagulation options, and vessel sealing, all within a compact, reliable generator.

This system also features autoclavable instruments, providing durability and extended use (up to 100 times).

We are excited to support your veterinary practice with these innovative solutions from Mano Medical via Genia.



CARBON TECH

CARBON TECHNOLOGY

- Energy is evenly distributed throughout the mat.
- Optimal temperature reached within a few minutes.
- Safer technology as it delivers a low contact temperature, while transferring enough energy to the patients to keep them warm. Does not store heat.


INTEGRATED TEMPERATURE PROBE

- Increased power for warming efficiently, while remaining safe thanks to the automated temperature regulation.
- Place the patient straight on the mat. Three speed levels: click on the power button to change the speed. Green & blue options are ideal for recoveries and red for surgeries.
- High sensitivity sensor and regulator: contact temperature maintained below 42°C (107.6°F).

HIGHLY RESISTANT MAT

- *Protective coating guaranteed for 5 years against bites and scratches.
- Easy cleaning and storage: power block can be disconnected.
- Safety shutdown after 12-hour-use.

YEARS
5
GUARANTEED*

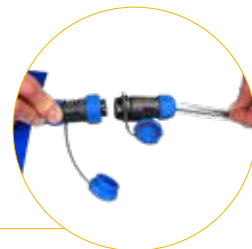

CARBON TECH

ADVANCED VET NURSE TRAINING, WA

Samantha, Owner

Did you know that rectal temperatures are not always the most accurate way to assess an anaesthetised patients core body temperature?

The patients rectum can dilate during anaesthesia, and the temperature of the rectal zone is not regulated as a priority by the animals body.

A better option for monitoring is the oesophageal probe, hands up who only just learnt that today like me! 🙌

I've recently been looking into heating options for patients under anaesthesia, and the safety concerns surrounding them and ive come across this carbon tech heating mat from Genia Australia.

It's designed to be quite soft and have no rigid heating elements so reduces the risk of compression sores and ischemia which can be an issue with older patients and those under anaesthetic for prolonged periods of time.

The other great aspect I love is the temperature probes to prevent it from overheating, even if the mat is used incorrectly.

Have you ever seen a burn from one of the old school style mats?

MT GRAVATT VETERINARY CLINIC, QLD

Renee, Vet Nurse

I have absolutely been loving the new carbon tech heat matts! Both sizes have been working really well in keeping all our patients warm under anaesthetic aswell as keeping them comfortable after procedures :)

The waterproof material and fast warming has been perfect for our daily surgery and hospital patients.

VETS ON PARKER, VIC

Kathryn, Vet Nurse

The heat mat is working really well so far, we use it for our sterile surgeries and we've found it's keeping our patients warmer for longer. The size fits perfectly to our table and it's easily to clean.

PYRAMID VET SURGERY, QLD

Siobhan, Vet Nurse

We love our Genia's heat mat, it heats up nicely and it's so easy to clean. We use it for all our surgeries.

ADELAIDE VET GOODWOOD ROAD, SA

Mandi, Nurse

Heats up quickly but does not overheat. Great that the power level can be adjusted.

Code	Wattage	Size	Description
800000	25W	52 x 62 cm	CARBONTECH MAT
800002	50W	52 x 112 cm	CARBONTECH MAT

SMART ✓ SAFE



Instant and even **heating**



Automatic power cut out



3 heating levels
50, 80 & 100%



Detachable cable
for easy cleaning

Increased safety with integrated regulation device


CarbonTech is equipped with a temperature probe and a high accuracy limiter, to ensure that the contact temperature does not exceed the safety limit (42°C).

Extra flat, Carbon Tech V2 carbon fabric does not store heat.

For higher safety, the power is turned off automatically after 12 hours of use.



42°C (107.6 F)
maximum

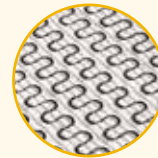

High density
carbon fabric

Carbon Tech technology

Carbon Tech incorporates a high density carbon fabric, specially designed for even heat distribution.

The soft and resistant coating offers more comfort to the patient.

Long-life and high efficiency

CarbonTech reaches its regulation temperature within just five minutes.

CarbonTech is designed to maintain the temperature necessary for warming the animal during long periods.



Strong and hygienic

A fastening system allows you to connect the power block safely and securely to the heat mat. The power block can be stored under the surgery table to ensure a safe working environment.

The protective coating of CarbonTech is designed for maximum resistance to animal scratching and biting and is totally waterproof.

CarbonTech can be easily cleaned with your normal veterinary cleaner. You can disconnect the power cable when cleaning.



Easy cleaning
with disconnectable
power block



WARMING SYSTEM

INSTRUCTIONS FOR USE:



- 1 Place your warming system on the surface where it will be used, such as a surgery table, recovery cage, etc. Ensure the mat is flat and positioned correctly to provide optimal warmth to the animal.
- 2 Connect your mat to the power block. To do this, align the two white dots on the plug. Once aligned, clip the plug together. Then, twist the blue part around the connection to secure it firmly. This ensures a stable connection that won't come loose during use.
- 3 Plug in the cable. Connect the cable from your warming mat to a power outlet. The fastening system allows you to securely attach the power block of the mat beneath the surgery table or around the cages. It provides a safe and efficient work environment without cluttering your workspace.
- 4 To turn on your mat, press the power button for 2-3 seconds. The mat will start warming up. Allow around 5 minutes to become fully operational.
- 5 Once the mat is on, you can select the desired power level. Press the same button once to cycle through the power levels. The color of the indicator light will show the selected power level: 100% (red), 80% (blue), or 50% (green). These power levels correspond to the amount of heat supplied by the mat, ensuring the animal receives the appropriate amount of warmth.

By following these steps, you can ensure a safe and effective warming process for your veterinary procedures, enhancing the care and comfort of the animals under your supervision.

EQUIPMENT



CARBON TECH

ABERDEEN VET CLINIC, VIC

Sharron, Practice Manager

We purchased the Genia warming System about 12 months ago for our wet / dental table. It's been wonderful, we have found our patients are staying warmer while under for long procedures and therefore recovering a lot quicker and happier.

RIPLEY VET HOSPITAL, QLD

Tanya Kunzel, Practice Manager

Easy and efficient to use with patient directly on the mat. Safe to use for our patients and easy to keep clean. Works amazingly at maintaining temperature during surgery.

TESTIMONIALS

CARBON TECH

BOORGAOON VET CLINIC, WA

Bree, Head Vet Nurse

We have been using our large carbon tech heating mat for a few months now and love it! It's super easy to clean, heats up quickly and does a great job at maintaining a patient's temperature during surgery. The padding is thick enough to support the patient, but not too bulky, so it's still easy to store.

KARDINIA VET CLINIC, VIC

Georgia, Practice Manager

All the staff really like the heat mat! It warms the patients really well and is easy to clean and store. Perfect for dentals as it just wipes clean!

TESTIMONIALS



VET DOPPLER

Blood pressure monitoring especially designed for pets

- Ergonomic design, small and light-weight
- Reliable, accurate, robust and liquid proof
- Rechargeable battery operated
- 4 cuffs in different sizes
- 3 different probes



♥ NON INVASIVE METHOD

Reliable & highly accurate

♥ FULL SET WITH A 8 MHz PROBE

with a large beam for better flow direction

♥ OPTIONAL PROBES

Pencil and transoesophageal

Adapted for multi-use



BLOOD PRESSURE PROBE

Large & sensitive beam for a quick flow detection



TRANSESOPHAGEAL DOPPLER PROBE

A cardiac monitoring solution during surgeries, on anaesthetised patients



PENCIL PROBE

For easy exotic animals' examination



DOPPLER SET INCLUDES:

- ✓ 1 8Mhz probe
- ✓ Battery charger
- ✓ Sphygmomanometer
- ✓ 4 cuffs in different sizes (4-8 cm, 6-11 cm, 7-13 cm and 8-15 cm)
- ✓ Headset
- ✓ Tube of gel

All the accessories are packed in a carrying case.

Code	Description
801000	DOPPLER BP - Full kit

SPARE PARTS

Code	Description
801012	Plastic cover
801002	Headset
801004	Pencil probe (exotics)
801006	Blood pressure probe 8mhz
801008	Oesophageal probe for monitoring
801010	AC adaptor for AUS
801014	Battery
801016	Cuff 4-8cm
801018	Cuff 6-11cm
801020	Cuff 7-13cm
801022	Cuff 8-15cm
801024	Sphygmomanometer



DOPPLER

BAYSIDE VETS, NSW

Naomi Illguth, Hospital Manager

The doppler is GREAT, works just as we hope. Nice and straightforward to use and clear sounds, we have no issues with it at all :)

VET PRAC, QLD

Margie, Owner

Genia Doppler has the same setup and application process as the Parks Doppler. The slim shape and size means that it is easy and convenient to use in the clinical setting. I noted that the probe was quite streamlined and the crystal is covered by a sturdy protective sleeve.

Lesley, Nurse

We use the VetTech BP doppler in our practice. We are happy with its performance, particularly because we can hear the patient's pulse throughout the procedure.

The improvements that we could recommend would be that the velcro on the cuffs be lengthened to keep the cuff in place more effectively.

BONDI VET HOSPITAL, NSW

Brooke, Nurse

This doppler is honestly the best and easiest doppler I've ever used. Placement and securing takes less than two minutes, even on the smallest of patients or the weakest of pulses. It's a nice small size, lightweight and the sound is super clear.

Keeli, Nurse

The DopplerVet is by far the most efficient and effective doppler I have used, with its compact storage case it is easy to ensure that when it is not being used it is protected. The battery life is long lasting and allows it to be in use the majority of the day without needing to be put on the charger in between patients. The probe is very efficient with finding the pulse and the sound is not dramatically disturbed by the smallest of movements.

Elsa, Nurse

The doppler has been great to use. The portable size allows easy transfer with patients. The wide beam probe makes it easier to locate the pulse and there is minimal background noise.

ANIMAL EMERGENCY SERVICE UNDERWOOD, QLD

Dr Jaxon Anderson, Vet

Nurses team I have been testing in clinic and it seems to be working exceptionally well for our critical patients. We tested it against the current doppler that we own and on a very sick patient with a very low blood pressure, the Genia was able to pick it up which is fantastic since the one we own currently could not pick it up.

TESTIMONIALS

DOPPLER

THE CAT CLINIC GOLD COAST, QLD

Dr Nicki, Vet

We've been using it for about a month now, I really like how it's compact, and feels quite sturdy compared to the Parks Doppler I'm used to. The cables, the crystal, and even the doppler device feels rock solid. The ability to plug in a normal headphone jack without added equipment is a plus too. Honestly I've got no complaints! I'm just hoping it doesn't break down as often as the Parkes and I'd be the happiest person alive!

ARH ESSENDON FIELDS, VIC

Danielle, Head Vet Nurse

The doppler is amazing. Everyone loves it. Picking up a pulse is super quick. The probe is robust, as expected.

ZOE MARKHAM, TAS

Zoe Markham, Training and consulting

I have used quite a few monitoring machines over my years of work. And although multi parameter monitors are fantastic, I never felt comfortable with the readings I got from cats. I much prefer Doppler reading for our feline friends. The Doppler Vet BP from Genia gives me everything I know and love in a Doppler. Clear headphones of you need to keep the room quieter (especially when taking a feeling blood pressure!), or an integrated speaker that had lovely coat done, multi sizes of cuffs for animals of all sizes, not just cats, the transversal linear probe has a lovely long cord so you don't have to be right next to the patient if you don't need to and its very sensitive probe allows for clear readings, the sphygmomanometer (or the squeeze BP display) is very sensitive to touch giving you great control over the pressure in the cuff and allows a quick release. I find this to be very helpful especially if your patient becomes uncomfortable for any reason and allows for multi readings (to get an average) to be done smoothly and efficiently. I also love that the Doppler comes with a thick gel which is imperative for getting readings. The machine itself is a breeze to use with straight forward smooth touch buttons. While Doppler's can be fiddly in the beginning to get used to the sensitivity and ease of use this machine makes using it a breeze. This is something you want to have in your monitoring kit not just for blood pressure but for checking blood flow in damaged limbs or during CPR when checking for blood flow for compression efficiency. It's a very handy unit to have in the clinic.

TESTIMONIALS

DOPPLER - PENCIL PROBE

EXOTIC VET NURSE TRAINING, WA

Samantha Clutterbuck, Owner

For reptiles, I've recently discovered these dopplers with pencil probes from Genia Australia that are amazing. The pencil probe can be used in place of a regular small square probe but I find it's especially helpful when looking for chelonian (turtle and tortoise) heart sounds.

The best place to find these heart sounds is generally the left thoracic inlet, but if you've ever tried this yourself, you'll know it can be quite deep, and require pushing your fingers onto the little square probe quite firmly to hear it.

TESTIMONIALS





Blood pressure measurement in veterinary medicine

At the beginning of the disease, the symptoms of arterial hypertension in animals are difficult to detect. In many cases, the diagnosis is made too late. Primarily older animals and particularly cats are affected by hypertension, however there are also a lot of indications seen in dogs.

A periodic check should be performed on older animals: blood analysis (detection of kidney insufficiency) and blood pressure measurement to detect hypertension. According to current studies, up to 30% of cats suffering from chronic kidney insufficiency have systemic hypertension. Cerebral, ocular, and severe renal problems can be a consequence of arterial hypertension.

It is also useful to perform blood pressure screening before or during an infusion: the arterial pressure can be a suggestive indicator of an animal's condition before an infusion. The veterinarian can adjust the infusion according to the results.

Reliable and highly accurate device

Measures blood pressure on small animals (including exotics) and detects blood pressure irregularities and changes.



Monitoring with the Doppler

The Doppler allows the flow detection in the artery. The probe can be affixed to the animal's paw or tail during an anaesthesia or during intensive care. The veterinary team will hear a sound signal reflecting heart rate and pulse.

INSTRUCTIONS FOR USE:

- 1 Select the appropriate cuff, attach it to the sphygmomanometer, and place it around the chosen artery (the radial, cranial tibial, and median coccygeal arteries are the most practical to monitor).
- 2 Clip the hair as short as possible or wet the hair with alcohol over the chosen artery.
- 3 Apply gel on the probe and place it in contact with the skin over the chosen artery. Don't apply too much pressure, as it could stop the blood flow. The part with the crystal should be placed parallel to the artery's blood flow.
- 4 Hold the probe softly with your finger or attach it to the limb with an adhesive bandage. The probe is correctly placed when you hear a clear pulsatile blood flow.
- 5 Inflate the cuff with the sphygmomanometer until the pulsatile blood flow is no longer audible.
- 6 Slowly deflate the cuff using the back trigger on the sphygmomanometer while observing the gauge.
- 7 Record your blood pressure reading at the first audible arterial pulse. This is your Systolic Arterial Pressure (SAP).
- 8 As the cuff deflates, listen carefully to obtain the Diastolic Arterial Pressure (DAP).
- 9 To validate the results, it's recommended to record multiple consecutive measurements, average the middle 3 readings, and record your average SAP.



CHOOSING THE RIGHT CUFF

The cuff should be approximately 30–40% of the circumference of the cuff site. The cuff size should be noted in the medical record for future reference.

It can be placed on either a limb or the tail, depending on the animal's anatomy and the clinician's preference. It is important to document the chosen cuff placement site in the medical record.

RECORDING MEASUREMENTS

Typically, the first measurement is discarded. Record at least 3, but preferably 5–7, consecutive and consistent measurements (with less than 20% variability in systolic values). Calculate the average of these values to determine the SAP. If there is any uncertainty, repeat the measurements.

Ensure you use the same limb and the same size cuff for all measurements. Record these details, including the multiple measurements, on your records form.



SEAL & CUT

Genia SEAL&CUT is designed to meet a wide range of surgical needs, offering both monopolar and bipolar options with various cutting and coagulation modes combined with vessel sealing.

- Time-out Function (Optional): automatically stops HF (high frequency) output after prolonged unintended activation.
- Toggle Foot Switch: allows switching between monopolar and bipolar modes using a dual-paddle footswitch.
- BiCoag Alarms: audio feedback alarm that alerts you when bipolar coagulation is complete, helping to reduce tissue charring and sticking to the forceps, to prevent overburning.
- Remote Power Setting: hand switch pencil can be used to control power remotely.



Code	Description
803000	MV1 full kit
803002	Rotation kit
803004	Autoclavable patient return electrode plate
803006	Autoclavable patient return electrode plate - cable
803008	Autoclavable bipolar forceps 20cm 1mm tip
803010	Autoclavable bipolar forceps 20cm 1mm tip - cable
803012	Autoclavable monopolar handswitching pencil
803014	Autoclavable sealing clamp 19cm 3mm tip
803016	Autoclavable sealing clamp - cable
803018	Autoclavable monopolar electrode – needle 5cm
803020	Autoclavable monopolar electrode – ball 5cm
803022	Autoclavable monopolar electrode – angled blade 5cm
803024	Autoclavable monopolar electrode – straight blade 5cm
803026	Autoclavable monopolar electrode – 8mm diameter loop 5cm
803028	Laparoscopic latching handle
803030	Autoclavable laparoscopic latching handle - cable
803032	Autoclavable maryland jaw sealer 18cm - insert
803034	Autoclavable straight fenestrated sealer 18cm - insert
803035	Autoclavable shaft 18cm
803040	Autoclavable laparoscopic sealer & divider kit - 18 cm

KEY FEATURES

Smart ESU - Electrosurgery unit

An advanced device that uses real-time tissue feedback to automatically adjust power output, ensuring precise cutting and coagulation.

- + Minimises tissue damage and surgical complications.
- + Enhances precision for safer, more effective surgeries.

6SENSE™ technology

An advanced feedback-controlled technology that maintains set power by sensing six parameters (voltage, current, power, tissue density, return electrode, and leakage RF current) at 4000 times per second.

- + Minimises electrode drag across tissues.
- + No need to adjust power setting as tissue changes.
- + Reduces thermal damage by 50%, promoting faster recovery.
- + Minimises collateral tissue damage and neuromuscular risks.
- + Decreases interference with other OT (operation theatre) equipment.
- + Reduces charring and sparking, ensuring cleaner cuts.

Smart ENDO-CUT mode

It divides cutting into cycles of cut and pause, automatically adjusting pauses based on tissue density. Optimised for polypectomy and papillotomy,

- + Enables precise cutting with coagulation in each pause.
- + Adapts cut-pause intervals to tissue impedance.

Randomised spray coag

Unlike conventional spray coagulation, this feature randomises spray patterns to achieve even coagulation across larger areas.

- + Sparks initiate further from tissue for uniform coagulation.
- + Lower power requirements, less necrosis, and faster results.
- + Enhanced performance when combined with Argon-enhanced ESUs.

BiCoag alarm

Audio feedback alarm that alerts you when bipolar coagulation is complete, especially useful in endoscopic procedures where tissue isn't directly visible.

- + Prevents over-coagulation and tissue sticking to the forceps.
- + Minimises smoke, tissue charring, necrosis, and prevent overburning.

SEAL & CUT

NICKLIN WAY VET SURGERY, QLD

Jasmin, Head Nurse

We're really enjoying using Genia Seal & Cut. We find the sealer especially helpful for spays, as it clamps easily onto fat. We've also tried the new cut/sealer handpiece during a splenectomy, and we loved it—it's already making a difference.



INCLUDED:

**SEAL&CUT - UNIT**

The main unit that features an all-in-one system for cutting, coagulation, and vessel sealing, complete with control panels and a touchscreen interface.

**MONOPOLAR FORCEPS**

Designed for precise cutting and coagulation, these forceps transfer electric current through the tissue, suitable for a range of surgical tasks. An option for a handpiece finger switch is available for cutting and coagulation modes.

**BIPOLAR FORCEPS**

These forceps offer greater control and safety by limiting the electric current between the two tips of the forceps, reducing the risk of injury to surrounding tissues during coagulation.

LOOP ELECTRODE

5cm, 8mm. Ideal for removing soft tissue or tumours, providing a precise cutting loop that minimises damage to surrounding areas.

ANGLED BLADE ELECTRODE

5cm. This electrode provides enhanced manoeuvrability for cutting in hard-to-reach areas, offering more precision in confined spaces.

BALL ELECTRODE

5cm. Designed for coagulation to control bleeding, it applies heat to a broader surface during surgery.

BLADE ELECTRODE

5cm. Perfect for making clean, straight incisions in tissues, commonly used in precise surgeries.

NEEDLE ELECTRODE

5cm. Excellent for pinpoint accuracy, this electrode is used in fine dissection or cauterisation of small, intricate areas.

**SEALING CLAMP**

Designed to effectively seal vessels and tissue, combining compression with energy to create a secure, lasting seal.

**LAPAROSCOPIC HANDPIECE**

The main control tool in laparoscopic surgery, designed to hold various inserts (like the Maryland jaw or fenestrated) for versatile use, offering enhanced precision and ease of operation during minimally invasive procedures.

FENESTRED INSERT

This insert guarantees a firm yet gentle grip on tissue, reducing trauma during laparoscopic procedures.

MARYLAND JAW INSERT

This insert is used in minimally invasive surgeries for delicate dissection and tissue handling, with a curved jaw for precise control.

**NEUTRAL PLATE (RETURN ELECTRODE)**

This component ensures patient safety by providing a grounding point for the electric current, helping to prevent burns by safely completing the electrical circuit during monopolar electrosurgery.

**FOOT SWITCH CONTROL**

Used to activate the electrosurgical unit hands-free, providing greater control and focus during surgery by allowing the surgeon to engage cutting or coagulation functions with a foot pedal.



MONOPOLAR ELECTROSURGERY



In monopolar electrosurgery, only one pole—the active electrode—is present at the surgical site. An electric current is delivered from the generator to the surgical site through the active electrode and returns to the generator via the patient return electrode. Monopolar electrosurgery is commonly used for general surgical procedures, with the return electrode (neutral plate) providing a safe path for the return of the RF current.

4 Monopolar cut modes

- **Low Cut** - used for cutting without sparking; useful in laparoscopic surgery for precise and clean cutting of thin tissue.
- **Pure** - the default monopolar cut mode, ideal for precise cutting without hemostasis.
- **Blend** - used when slower cutting and moderate hemostasis are desired.
- **Endo** - pulsed cutting mode with an appropriate pulse interval.

4 Monopolar Coag Modes

- **Soft (Desiccate)** - desiccation dehydrates and destroys tissues without sparking, as the active electrode is in direct contact with the tissue; ideal for coagulating soft and delicate tissues.
- **Swift** - soft coagulation covering a larger tissue area.
- **Fulgurate (Low)** - the default monopolar coag mode, suitable for non-contact coagulation with moderate sparking across various applications.
- **Randomised Spray (High)** - designed for distance coagulation over large areas with minimal tissue damage and necrosis.

BIPOLAR ELECTROSURGERY



In bipolar electrosurgery, both the active and return electrodes are located at the surgical site, so a separate patient return electrode is not required. The electric current flows from the active electrode to the return electrode through the tissue grasped by the instrument. This technique is mainly used in delicate surgeries, cosmetic procedures, and neurosurgery.

2 Bipolar cut modes:

- **Macro** - provides precise bipolar cutting with hemostasis.
- **Pulse BiCut** - allows clean cutting while maintaining hemostasis.

4 Bipolar coag modes:

- **Micro (Precise)** - maintains low voltage to avoid sparking; ideal for delicate bipolar tissue desiccation.
- **Standard** - general bipolar coagulation for various tissue types with moderate hemostasis.
- **Force** - designed for bipolar coagulation of thicker and high-impedance tissues.
- **Sealer** - effective for standard vessel sealing and fusion.



INSTRUCTIONS FOR USE:

Let's get you started with your Seal & Cut MV1. : your electrosurgery unit which allows for monopolar and bipolar electro surgery, combined with vessel sealing.

Ask
your
guide
today

- 1 Connect the foot pedal to the back of the unit. Make sure the small indentation on the plug faces up when connecting. To connect your monopolar pen, you will insert its very unique cable into its designated monopolar connection point. To connect your bipolar forceps, you will insert its cable to its exclusive bipolar connection point. To connect your vessel sealing clamp or your laparoscopic handpieces, you'll insert their specific cable to its unique vessel sealing connection point. The neutral plate connects into the connection point labelled "patient".
- 2 Plug in the device and locate the power button at the bottom left of your device. Press it, and the button will turn green. You will hear a welcome sound, followed by a screen with three patient size pre-set options: small dog/cat (<10 kg), medium dog (10-30 kg) , and large dog (>30 kg) Different skin types require different power levels, so the pre-set options provide a good starting point. Begin with the option corresponding to your patient, and adjust as needed. Once you select the animal size, a new screen will appear, showing four quadrants:
 - Monopolar cut
 - Monopolar coagulation (COAG)
 - Bipolar cut
 - Bipolar coagulation & vessel sealing

To adjust the power in any of these modes, simply tap on the screen to the corresponding quadrant and use the up and down arrows to set the desired level. After adjusting, press the back arrow to save the changes.
- 3 The foot pedal controls whether you are using monopolar or bipolar functions. You can see its position on the screen. If you're using monopolar, the foot pedal must be on the top line (on your screen) and for bipolar or vessel sealing, on the bottom line. Switch between modes by pressing the black button in the centre of the foot pedal.
- 4 The vessel sealing function is an extension of the bipolar coagulation mode. A seal is achieved by passing bipolar current through the tissue, causing the collagen and elastin to fuse together without carbonisation. This function appears in the bottom right quadrant of the screen, which will be blue when set in bipolar mode. To activate vessel sealing, select the quadrant and switch the mode to vessel sealing – the colour will turn purple, indicating that the vessel sealing mode is now active.

AREAS OF APPLICATION:

APPLICATIONS OF ESU / DIATHERMY SECTION

- Gynecology
- Urology (Resection)
- Laparoscopic Surgeries
- Endoscopy (ERCP)
- Oncosurgery
- Heart & Chest Surgery, CABG
- Neurosurgery, Spine
- Arthroscopy
- All Other General Surgical Procedures

APPLICATIONS OF VESSEL SEALING SYSTEM SECTION

- Lap-Assisted Vaginal Hysterectomy (LAVH)
- Vaginal Hysterectomy
- Gastrectomy
- Gastric Bypass
- Appendectomy
- Colectomy
- Nephrectomy
- Nissen Fundoplication
- Colon Resection
- Cystectomy
- Radical Hysterectomy
- Liver Resection
- Adrenalectomy
- Splenectomy
- Salpingo-Oophorectomy
- Radical Prostatectomy
- Abdominal Hysterectomy
- Hemorrhoidectomy



vétopop

The best for my pets



The BEST for my PETS !

Vétopop products are designed for veterinary clinics and pet owners who want to maximise the joy and enrichment they share with their furry companions. We know that dogs and cats thrive on a healthy balance of physical activity and mental stimulation.

That's why Vétopop offers a range of stimulating and fun products that promote movement, learning, and overall well-being.

Our products are designed to complement veterinary services and support key aspects of pet health:

- **Weight management:** slow feeder bowls to encourage mindful eating
- **Boredom relief:** Pop ball and sea urchin toys for interactive play
- **Dental health:** dental toys designed to support clean teeth and gums
- **Cooling vests or coats:** to keep pets comfortable and cool / warm during warmer / cooler months

Vétopop products provide great added value for both veterinary clinics and pet owners, helping to enhance pets' lives and support their overall wellness.

Products available through veterinary clinics only.

Dental Toy

Rotating rubber dental toy for dogs. Durable dog toothbrush improving oral health.

Textured bumps and ridges help removing dental plaque.

Can be disassembled to add some treats or take some rings off when cleaning.

Reduces boredom: the dental toy keeps your animal entertained.

Made of durable elastomer: light and sturdy.

Easy to clean & dishwasher safe.

Sold individually.



👍 Dishwasher safe !

Length 14 cm
Diameter 6 cm

● code 700450

TESTIMONIALS

DENTAL TOY

Jacqui & Arlo, Pet Owner

The dental chew toy is really good for her because she loves chewing and hasn't been able to destroy it yet which is saying something for her.



Treat dispenser puzzle

Wooden interactive puzzle feeder for cats and dogs.

Place dry food or treats in the compartments and let your pet play with the sliding blocks to reveal the food.

An interactive game that boosts pets' mental stimulation and keeps them engaged.

Can be used to slow down your pets eating pace.

Made of water-resistant composite wood fibre.

Handwash with soapy water.

Sold individually.



Round model
Diameter 25 cm

code 700510

Novelty model
Width 29 cm
Length 19 cm

code 700512



Crazy Urchin Toy

Treat dispenser for dogs and cats. Sea urchin shape gives unpredictable bounces.

Wide opening: easy to fill in with dogs' favourite treats.

Unpredictable bounces: encourage physical activity and prevent from boredom.

Bright, colourful & funky design to be easily found in the grass.

Made of durable elastomer: highly resistant.

Suitable for any breed: comes in 2 sizes, pink and yellow.

Easy to clean & dishwasher safe.

Sold individually.



👍 Dishwasher safe !

Size S (small dog/cat)

Width 5 cm

Length 6 cm

● code 700502

Size M (medium/large dog)

Width 7 cm

Length 9 cm

● code 700503



TESTIMONIALS

SEA URCHIN

MURRAY BRIDGE VET CLINIC, SA

Kirstie, Vet Nurse

Our border collie puppy Nova loves the high bounce and compact size of the urchin toy. Especially when there's food stuffed in it.



Bouncing ball

Treat dispenser bouncing ball for dogs and cats.

Keep your pet entertained for hours.

Prevents boredom and encourages physical activity.

Can be used with treats.

Made of durable elastomer.

Sold individually.



👍 Dishwasher safe !

Size S (small dog)

Diameter 6 cm

● code 700400

Size M (medium/large dog)

Diameter 9 cm

● code 700402

TESTIMONIALS

POP BALL

MURRAY BRIDGE VET CLINIC, SA

Kirstie, Vet Nurse

Our puppy loves the Vétopop Bouncing Ball. She's often napping on it or just carrying it around. It's the only ball she will actually go fetch (won't bring it back though!). Definitely one of her favourites and is still in perfect condition despite her needle sharp teeth (she's only 16 weeks old).

Pop lick pad

Square shape with cross pattern. Strong suction cups underneath to stick the pad on surfaces.

Licking improves pet's oral health, as it helps their digestion.

Stimulates animals mentally and prevents from boredom.

Distracts patients during consultation (e.g. puppy vaccinations) and grooming.

Can be placed in the microwave and the freezer. Can be used with cold and hot food.

Dishwasher-friendly.

Sold individually.

**👍 Ideal for puppy schools.
Can be used on your consultation tables.**

code 701400



TESTIMONIALS

LICK PAD

WINGHAM & VALLEY VETS, NSW

Amber, Head Nurse

We love the lick pads and use them in all our consults to promote a positive association with vet exams for our patients. The suction cup feature is great and we sell quite a lot to our clients as they find them great too !





TRAVEL BOWLS

Two collapsible bowls with carrying case. Keep your furry friends hydrated and happy on your travels.

Made of non-skid silicone.

Velcro layer underneath for easy removal.

Dishwasher safe.

Space-saving: easy to carry and store in your car, luggage, or backpack. Ideal for travel.

2 bowls are better than 1! No need to choose between water and food, pets can have both simultaneously.

Sold individually.



**Easy to clean,
dishwasher-safe!**

750 mL

● code 700304



Slow feeder bowls



**Slow down eating
pace up to 10 times!**




Improve digestion



Boost mental stimulation



**Ideal for overweight
pets and overeaters**

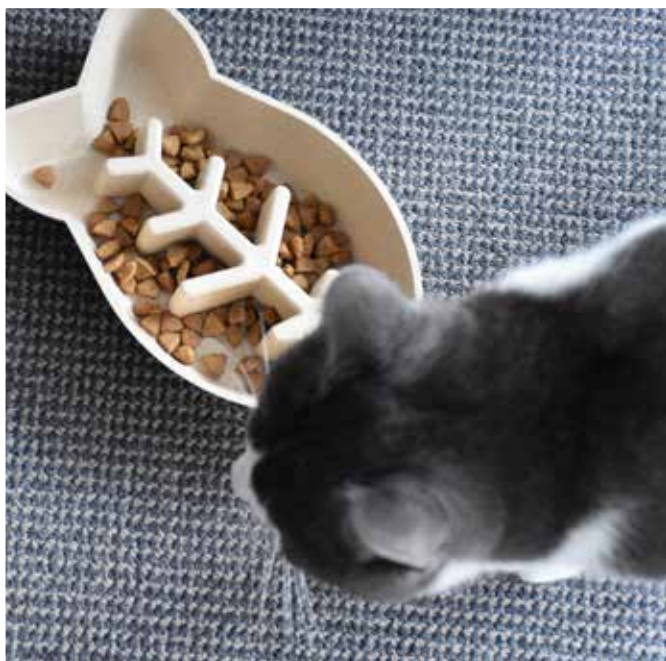

BIO Slow feeder bowl

The eco-friendly slow feeder bowls for cats and dogs.

Plant-based material: made of vegetal fiber that comes from sustainable forests.

Helps reduce greenhouse gas emissions.

100% biodegradable.

100% compostable.

Handwash only.

Packaging made of recycled & recyclable cardboard.

Sold individually.



Eco-friendly.
Made of plant-based plastic.
100% biodegradable.
100% compostable.



Bone-shaped, size M
(medium/large dog)

1200 mL

code 700312 ●

Bone-shaped, size S
(small dog)

840 mL

code 700310 ●

Fish-shaped
(cat)

350 mL

code 700308 ●

Slow feeder bowl

Slow feeder bowls for cats and dogs. Keep your pets engaged during meal time and slows down their eating pace. Reduce gulping, for a healthier eating behaviour. They improve digestion and reduce bloating.

Slow down eating pace up to 10 times: ideal for overweight and bloated pets.

Boost mental stimulation.

Gripping pads underneath.

Dishwasher safe.

Sold individually.



Easy to clean,
dishwasher-safe!



Bone-shaped, size M
(medium/large dog)

1200 mL

code 700302 ●

Bone-shaped, size S
(small dog)


840 mL

code 700299 ●
code 700300 ●
code 700301 ●

Fish-shaped
(cat)


350 mL

code 700305 ●
code 700306 ●
code 700307 ●

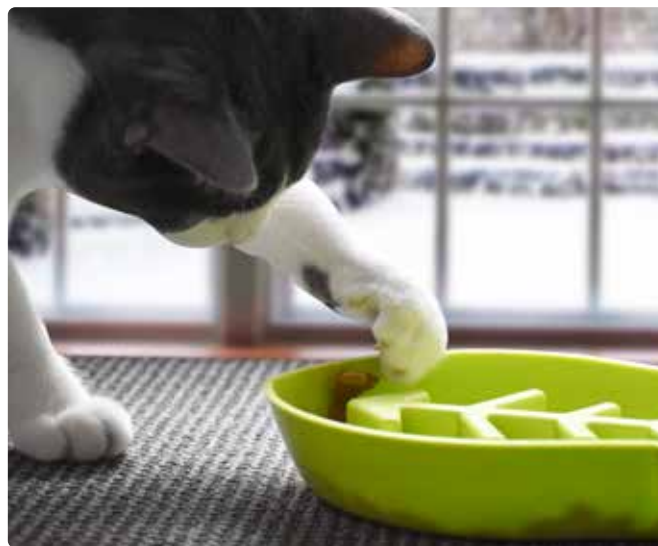



Cosypop


The ideal surgical and recovery bed!

 **MULTI-PURPOSE** – Ideal for SURGERIES: it provides insulation and protects patients from cold surfaces. The COSYPOP® is particularly helpful for airway procedures, dental surgeries and orthopaedics, thanks to its thick non-absorbent headrest. Ideal for in-clinic RECOVERIES: the bed reduces pressure on muscles, joints, and prevents from bedsores. It is non-absorbent: liquids go through, keeping animals dry and clean.

Can also be used as a PET BED at home or as a PUPPY BASKET.

 **MATERIAL** – Unique meshed material that makes the bed breathable and non absorbent. It provides insulation and keeps animals dry at all times. The headrest and the base are made of the same components.

 **NON-SLIP BASE** – Anti-slip base to keep the COSYPOP® stable on surgery table, in the cage or kennel.

 **MACHINE WASHABLE** – The base and the headrest can be washed separately by unzipping them. Up to 40°C on gentle cycle. Quick drying. Do not tumble dry.

IDEAL FOR YOUR:

In hospital recoveries

Dental procedures

Puppies basket

Surgery tables

Pet bed



40°C - Do not tumble dry
Wash in delicate cycle
if using top loaders



TESTIMONIALS

COSYPOP

BRISBANE VETERINARY SPECIALIST CENTRE, QLD

Dr Kate, Vet

The Cosypop has been a hit with the team. It's extremely useful for patients with poor bladder control, to help drain away moisture and avoid urine scalding. The soft cushioning in the base is also very comfortable for our non-ambulatory patients, and the raised edges allow us to easily keep their head elevated. It's super simple to clean, and fits perfectly in a standard size cage. We often reach for it if we have patients with tick paralysis, a neurological disorder, or just any animal that needs a particularly comfy spot that feels like a more homely dog bed.






Machine washable !



2 IN 1 !




Size 1

33 x 50 cm

code 701300 
code 701302 
code 701304 

Size 2

39 x 71 cm

code 701306 
code 701308 
code 701310 



Cooling vest

Dog cooling vests are a great way to ensure the safety of dogs during hot, humid conditions. When drinking cold water and shade are not enough, a cooling vest can provide extra comfort to ward off heat exhaustion or heatstroke.

Dunk the vest under cold water, wring it out and put it on the dog: the micro-perforated cooling fibres will keep the dog cool for several hours.

Improves animal well-being and reduces stress.

Do not leave the vest on the dog once it is dry.

Velcro fastening.

Reusable.

Hand wash only.

Sold individually.



Size S	Size M	Size L	Size XL
A Length 35 cm	Length 40 cm	Length 50 cm	Length 60 cm

● code 700700 ● code 700701 ● code 700702 ● code 700703



COOLING VEST

LATROBE VET GROUP, VIC

Sally, Operations Manager

I used it for my bulldog last year, it's fantastic! It helped him to keep cool all along the summer.

TORQUAY ANIMAL HOUSE, VIC

Ann, Practice Manager

We use the cooling vest for overheating patients. It works perfectly as it provides the right cooling amount for the animal.

CANNING VALE VET, WA

Sophie, Nurse

This dog cooling vest is awesome and works extremely well. Unlike the cooling collars that don't fit around my pugs necks, the cooling vest fits around their bodies perfectly. I put it in water, wring it out and they stay cool. I highly recommend it for all pets with this summer heat ☀️

CLAYFIELD VETERINARY CLINIC, QLD

Megan McCubbin, Nurse

My dog adores her cooling coat. She is normally not one to enjoy wearing coats in winter, however she finds it very beneficial when dealing with the harsh humidity of the Queensland climate. It is especially beatifical during an afternoon walk, or after a walk. She seems to settle down from the heat twice as fast and she look a lot more comfortable. I normally keep it on her for an hour and a half after doing something strenuous.



Reversible dog puffer

Reversible orange and grey puffer for dogs.

Thanks to its synthetic fabric padding, it protects your dog from cold temperatures.

Easy to put on thanks to its Velcro ventral fastening.


Perfect for safe walks at night thanks to its reflective strips around the zip.

Delivered with its own storage pouch to save space. Easy to carry.




Machine washable up to 30 °C.

Sold individually.

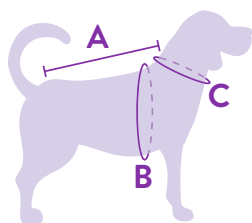 **Light & warm!**



Size XS	Size S	Size M	Size L
A : 27-30 cm	A : 38-40 cm	A : 42-45 cm	A : 44-47 cm
B : 43-45 cm	B : 52-54 cm	B : 62-64 cm	B : 69-72 cm
C : 30-33 cm	C : 33-38 cm	C : 35-38 cm	C : 41-44 cm

 code 700710  code 700712  code 700714  code 700716



NEW

TRAVEL POP

Foldable and practical carrier!

Foldable carrier made of rigid polypropylene with two large open side gates on the front for easy reach access. Ventilation along all sides.

5s only to assemble or fold - Easy storage.

Locking clips to keep your pet safe.

Top handle to carry your pet easily.

Easy to clean.

Recommended for pets (cats, small dogs, rabbits) up to 5 kg.

Sold individually.

 **FOLDABLE
& SAFE**

 **5s only to
assemble or fold**

 code 701500



VÉTOPOP



Dog poop bags dispenser

Silicone poop bags dispenser.

Refillable.

Can be easily attached to the leash.

Includes one roll of 20 disposable bags.

Sold individually.



● code 700414

Dog poop bags

Funky dog poop bags.

Rolls of 20 printed poop bags (10µm).

Made of polyethylene.

Sold in packs of 4 rolls of 20 bags.



● code 700800

TESTIMONIALS

POP BAG DISPENSER

EMERGENCY VETS 247, QLD

Dr Jacqueline, Vet

The poo bag dispenser is very ergonomic and portable - I like that it is soft rather than rigid, and also enjoy the dispensing design.



DRYPOP

Reusable training pad.

Reusable: Environmentally friendly.

Leakproof & highly absorbent.

Odour controlling.

Non slip layer.



Ideal for puppy training or in cages.

Sold individually.

 **ECOFRIENDLY, WASHABLE,
NON-SLIP**

Size 1	Size 2	Size 3
21 x 35 cm	45 x 60 cm	60 x 85 cm
80 mL	300 mL	600 mL

 code 700816  code 700818  code 700820

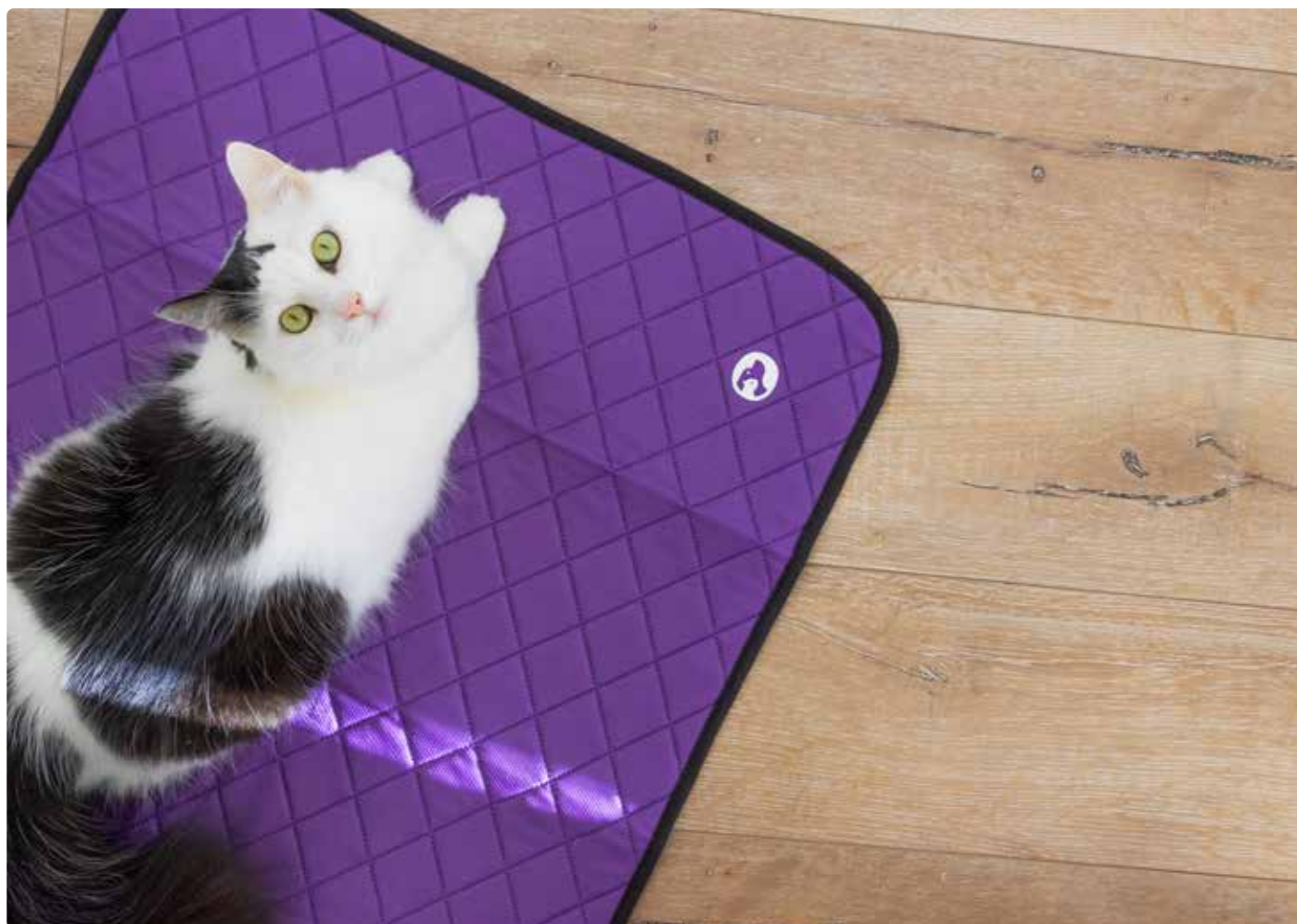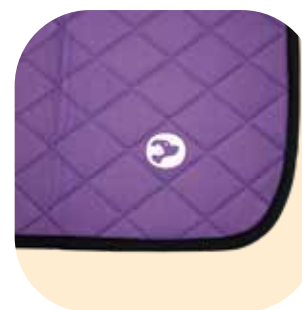





LARGE ANIMALS & EQUINE



PROTECTION

ECOGAN

Polyethylene GREEN gloves without neck strap (25 µm thickness), 95 cm length. These gloves are strong and thick providing good value exploration gloves, for day to day use.

Sold in boxes of 100 gloves.



RESISTANCE
2/4

PROTECTION
2/4

SENSITIVITY
2/4

Code x 100	Length
---------------	--------

307650 ●	95 cm
----------	-------

EXPLOGAN®

Polyethylene ORANGE gloves with a neck strap to prevent the glove from slipping (30 µm thickness), 95 cm length. Gloves known for their excellent puncture resistance with high sensitivity. Ideal for equine vets, requiring gloves with fine tactility and high penetration resistance.

Sold in boxes of 100 gloves.



IDEAL FOR
EQUINE VETS



RESISTANCE
3/4

PROTECTION
3/4

SENSITIVITY
3/4

Code x 100	Length
---------------	--------

307611 ●	95 cm
----------	-------

POLYGAN

Polyethylene ORANGE gloves with large armpit protection area (24 µm thickness), 124 cm length. Stretchable glove to prevent from slipping.

Sold in boxes of 100 gloves.



RESISTANCE
4/4

PROTECTION
3/4

SENSITIVITY
3/4

Code x 100	Size	Length
---------------	------	--------

307634 ●	S	87 cm
307636 ●	M	90 cm
307638 ●	L	95 cm

PROGAN®

Polyethylene YELLOW gloves without neck strap (30 µm thickness), 90 cm length. The tear resistance and sensitivity is optimal for obstetrics, rectal examination, as well as artificial insemination.

Sold in boxes of 100 gloves.



RESISTANCE
4/4

PROTECTION
3/4

SENSITIVITY
3/4

Code x 100	Size	Length
---------------	------	--------

307626 ●	S	87 cm
307627 ●	M	90 cm

TESTIMONIALS

PROGAN

SMITHTON VET SERVICES, TAS

Dr Craig, Vet

Genia exploration gloves are thinner than others brands, providing a better sensitivity. They don't tear apart and last long.



SENSIGAN®

Polyethylene CLEAR gloves with a neck strap to prevent the glove from slipping (28 µm thickness), 95 cm length. Gloves known for their excellent puncture resistance with the highest sensitivity, makes them ideal for mare reproduction work.

Sold in boxes of 100 gloves.



RESISTANCE

3/4

PROTECTION

3/4

SENSITIVITY

4/4

Code
x 100 Length

307610 95 cm

MAXOGAN

Polyethylene GREEN gloves with large armpit protection area (25 µm thickness), 124 cm length. The tear resistance is optimal and the shoulder protection make them the perfect gloves for reproduction work and rectal examinations.

Sold in boxes of 50 gloves.



RESISTANCE

2/4

PROTECTION

4/4

SENSITIVITY

2/4

Code
x 50 Length

307623 120 cm

POLYBOOT

Disposable plastic overboots with elastic bands. Ideal for farm visits.

Sold in packs of 50 boots (25 pairs).



Code Pack.

307640 Bag of 50



GÉNIAPLEX

Hypodermic plexiglas reusable syringes. The Géniaplex has a lifetime guarantee!

- Manufactured from a solid block of plexiglas, treated to avoid cracking.
- Stainless steel plunger shaft and washer: won't distort under any sterilisation.
- Metal top and bottom caps making it sturdier for heavy duty use.
- Resistant to boiling water and sterilisation.
- Can be entirely dismantled.
- Black indelible graduations.
- Bibloc luer-lock hub.

Sold individually with instructions for use, an extra washer and silicone oil.

HYPODERMIC SYRINGE

	Code	Model	Capacity	Hub	Graduation	Washer
A	100010	Standard	20 ml	Eccentric	1 ml	20-28
	100012	Short	30 ml	Eccentric	2 ml	24-32

DRENCHING SYRINGE

	Code	Model	Capacity	Hub	Graduation	Washer
B	180052	Drenching syringe	300 mL	Concentric	10 mL	40-50



IDEAL FOR
REPETITIVE USE



ELPLEX

Moulded plastic hypodermic syringes with engraved graduations.

- Easy to use, fit comfortably in your hand.
- Plastic plunger, cap and handle.
- Resistant to boiling water and autoclave up to 134°C.
- Broad finger grips for maximum control.
- Flat, beveled and extra-strong silicone washers: optimal glide of the plunger.
- Nickel-plated brass hub.

Sold individually.



	Code	Capacity	Hub	Graduation	Washer
	103002	5 mL	Concentric	0,1 mL	9-14
	103003	10 mL	Concentric	0,5 mL	14-20
	103504	20 mL	Eccentric	1 mL	20-28
	103505	30 mL	Eccentric	2 mL	24-32
	103506	50 mL	Eccentric	5 mL	30-40

SILIGLISS SILICONE OIL

Thick 10mL silicone oil to lubricate syringes and drenching guns washers.

Sold individually.



	Code	Type	Capacity
	106023	Tube	10 mL






IMOPLAST

Premium IV set with double needles for large animals.

- 150 cm clear PVC tubing.
- Plastic Luer-Record hub.

Sold in non-sterile bags of 25 IV sets.

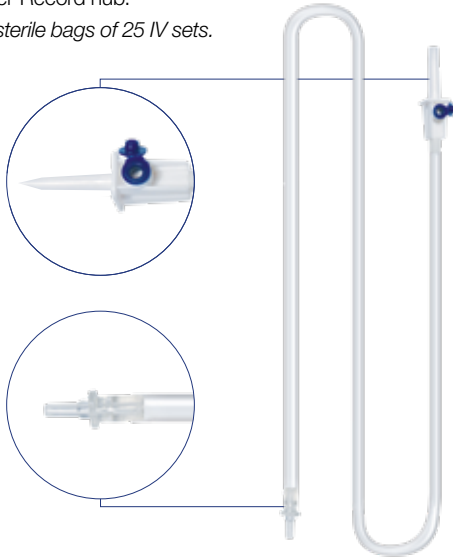

Code x 25	Description
140514	Complete with plastic Luer-Record hub

RUMEN NEEDLE

Stainless steel intra-rumen needle. Special needle for direct injection through the hide into the rumen of cattle.

- Fit a 4 mm tubing or connect it to any luer syringe.
- Tri-beveled, rounded heel.
- 160 mm x 4 mm - 8G x 6" 3/8.
- 6% taper, Luer-olive hub.

Sold in boxes of 12 rumen needles.




Code x 12	Length	Size (G)
147090	160 x 4 mm	8G x 6 3/8"

STAINLESS STEEL UTERINE CATHETER

Stainless steel uterine catheters for large animals.

- 540 mm long.
- Blunt, closed end with two side holes.
- Luer-Lock hubs.

Sold individually.

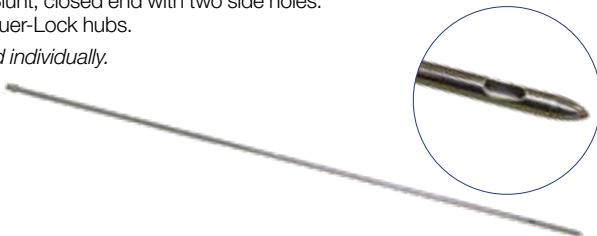


Code	Type	Length	Ext. diameter
307110	Straight (large animals)	540 mm	5 mm



CAESAREAN FORCEPS

Stainless steel caesarean forceps with serrated jaws for holding the uterus up to the abdomen opening. Replaceable red silicone jaws available.

Sold individually.

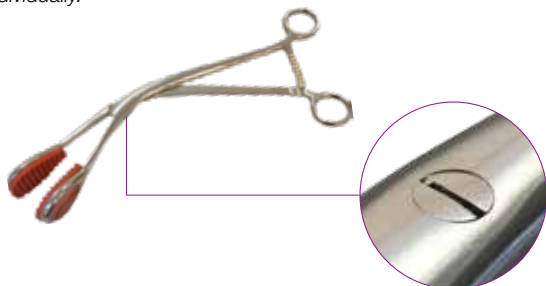


Code	Description	Length
227077	Screw	28 cm
227087	Replacement silicone jaws, pair	-

C-CUT

Single use plastic caesarean scalpel for large animals. Internal blade behind an extended plastic tip. Non sterile.

Sold in packs of 10 c-cut.

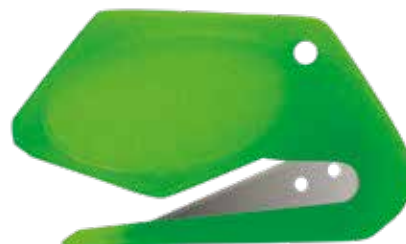

Code x 10	Description
307985	C-cut

BUHNER'S VAGINAL TAPE

Braided silk suture tape for preventing prolapse. Less dangerous and more secure than safety pins.

Sold individually.




Code	Length	Int. diameter	Use
220850	5 m x 7 mm	6,5 mm	Cow
220852	5 m x 4 mm	2,5 mm	Sheep

BUHNER'S NEEDLES

Stainless steel suture needles with handles, used to repair vaginal and uterine prolapses.

Sold individually.

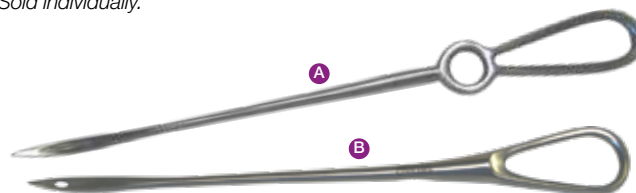

	Code	Model	Length
A	227900	Standard	30 cm
	227901	Tapered	30 cm
	227902	Standard	20,5 cm
B	227903	Tapered	20,5 cm




RUMEN TROCAR

Stainless steel rumen trocar.

- High quality trocar with plastic cannula and safety cap for emergency bloat.
- Trocar shaft 18 cm long x 9 mm diameter.
- Cannula with two holes.
- Extra cannula available.

Sold individually.

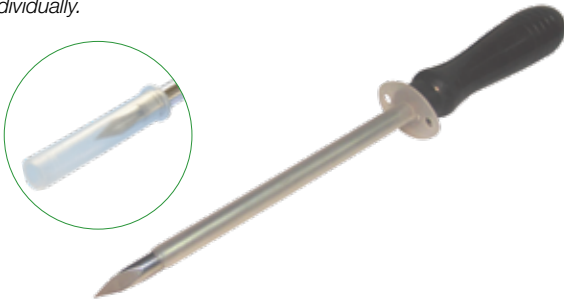

Code	Description	Int. diameter	Ext. diameter	Cannula length
220340	Rumen trocar with cannula	9 mm	10 mm	162 mm
220460	Cannula only	9 mm	10 mm	162 mm

SCREW TROCAR

Plastic screw trocar. Stainless steel pin ; cannula 5 mm.

- Helps animals suffering from chronic bloat. Especially useful for young stock. The cannula can be left inserted for up to three weeks, without arm (stitch it on).
- The threaded portion holds the rumen against the peritoneum, preventing peritonitis.

Sold individually.




Code	Description
220600	Complete trocar with instructions

SAW WIRE

Dehorning and embryotomy wire. 2 options: standard with 18 strands and swedish with 27 strands (sharper).

Sold individually.

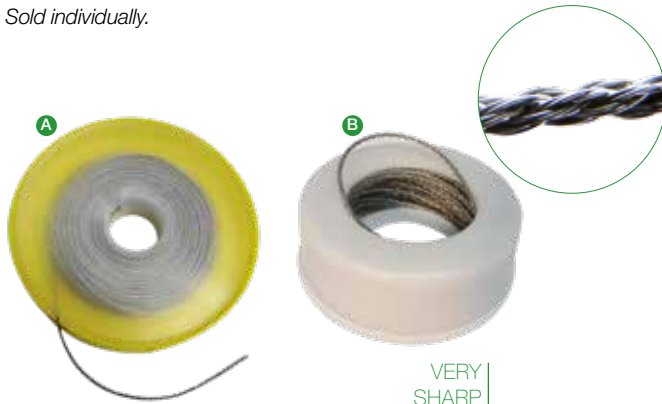

	Code	Description	Length
A	220840	Swedish 27 strands	12 m
B	227842	Standard 18 strands	12 m

SAW WIRE HANDLE

Saw wire handle. The knurled ends help to maintain a secure grip. Length: 10 cm.

Sold in pairs.



Code	Description	Length	Diameter
227840	Thin in stainless steel	103 mm	12/15

TESTIMONIALS

SAW WIRE

TAREE VET HOSPITAL, NSW

April, Practice Manager

Genia saw wire is good and reliable. It's also good value for money compared to others brands.

DANISH TEAT SLITTER

Stainless steel. Blades open by loosening the small burr on the side of the instrument. Open teat sphincters for better milk flow. Useful for follow-up surgery and mastitis.

Sold individually.



Code	Length
307984	19,5 cm

ECOMAG

High quality green double magnet to help prevent Hardware Disease in cattle. Multiple retaining bars.

Sold in boxes of 12 ecomags.



Code x 12	Description
227820 ●	Green



PENMARK®

High quality marker pens for cattle, sheep, goats and pigs. Marker pen holder frees up both hands.

- Velcro band for hand fitting.
- Machine washable Neoprene.

Sold in boxes of 10.

Code	Model
227628	Marker pen holder



PISTOLMAN

Automatic drenching guns with interchangeable barrels.

- Accurate and quick adjustment via the rear part.
- High grade and heavy-duty plastic material which can be entirely dismantled.
- Transparent and interchangeable barrels, resistant to boiling water sterilisation, with indelible & engraved graduations.
- Fitted with brass or chrome plated nozzle.
- Silicone o-rings.

Sold with a nozzle, instructions for use and silicone oil.



Code	Capacity	Graduation	Model	Nozzle model	Nozzle diameter
180300	20 ml	1 ml	Chromed nozzle	Short	6 mm
180065	50 mL	2 mL	Brass nozzle	Short	10 mm
180070	70 mL	5 mL	Brass nozzle	Short	10 mm
180075	150 mL	10 mL	Brass nozzle	Long	10 mm
180079	250 mL	25 mL	Brass nozzle	Long	12 mm
180078	250 mL	25 mL	Brass nozzle	Long	18 mm (only for thick products)

SYRMAN

Translucent plastic drenching syringe with indelible engraved graduations. Can be used as a horse mouth flusher.

- Use as a one-shot drench gun for bloat, glycerine, obstetric lubricant.
- Resistant to boiling water sterilisation.

Sold with a 19 cm long x 10 mm diameter nozzle and silicone oil.



Code	Capacity	Model
180050	300 mL	Brass nozzle
180057	450 mL	Brass nozzle
180336	450 mL	Chromed nozzle

FEEDING SYRINGES

Disposable 60 mL feeding syringes to administer emergency colostrum or electrolytes directly to the stomach.

- 1 mL graduations
- 31cm long x 6mm diameter feeding tube.

Sold in a kit including the feeding syringe and tube.

Code
107031



ORALPLEX

Drenching syringe.

- Precise dosage regulation with transparent barrel to facilitate viewing of liquids.
- Sturdy, non-distorting, heavy duty plastic material.
- Silicone washer to prevent leakage.
- 6 mm curved brass canula with rounded nozzle for safety (screw on).
- Resistant to boiling water and autoclaving.

Sold with a nozzle, instructions for use and silicone oil.



Code	Capacity	Assembly	Model	Nozzle diameter
180022	50 mL	Screw	Elplex syringe	6 mm





STOMAFLEX & NASOFLEX

STOMAFLEX **A**

Red silicone nasoesophageal tubes for horses and cattle.


- Silicone tube with open end and a side hole.
- Flexible and atraumatic.

NASOFLEX **B**

PVC nasoesophageal tubes for horses.

- PVC tube with open end and two side holes.
- One end widened for pump or funnel.
- Atraumatic and optimal flexibility for an easy insertion - Quick to clean.

Sold individually.

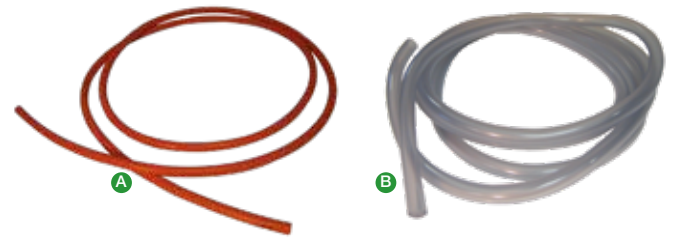

STOMAFLEX				NASOFLEX		
Code	Diameter ext.	Diameter int.	Length	Code	Ext. diameter	Length
227013	9,6 mm	6 mm	2 m	222000	9,5 mm	2,6 m
227005	12,6 mm	8 mm	2 m	222002	13 mm	2,7 m
227010	15,6 mm	9,5 mm	2,25 m	222004	16 mm	3,4 m
227014	19 mm	11,6 mm	3 m	222006	19 mm	3,6 m

NASOFLEX PUMP

Sturdy plastic pump with tube adapter for stomach tubes. Flow 0,5L.

Sold individually.



Code	Description
222020	Plastic pump + tube adapter

ESOFLEX

Red plastic-coated metal stomach tubes for adult cattle (18 mm diameter) and calves (14 mm diameter).

- Flexible rod with a ball at the end.
- Tool used to remove foreign bodies by reaching obstruction lodged in the oesophagus and gently pushing them into the stomach.

Sold individually.



Code	Type	Length	Diameter
227015	Calf	1,35 m	14 mm
227016	Adult	1,60 m	18 mm

DRENCHING NOZZLE

Screw model - brass drenching nozzle for hypodermic syringes.

- 12 cm length, 6 mm diameter.
- Olive connector.

Sold individually.




Code	Description
187036	Model with screw in place of the hub

DRENCHING WASHERS

Drenching washers for Génia drenching guns.

- High duty silicone washers.
- Water-resistant.
- Guarantees a perfect seal and a smooth glide of the plunger.

Sold in packs of 6 drenching washers.

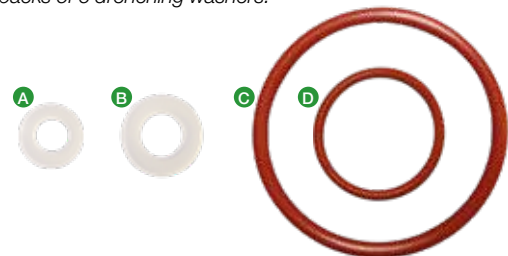

	Code	Description	Pack.
A	109267	Géniaplex/Géniamatic 20, 30 mL long	x 6
B	109268	Géniaplex 30 mL short, 50/60 mL long	x 6
D	189501	Syrman 300/450 mL	x 6
C	189514	Géniaplex/Géniamatic drenching gun	x 6



SURGIFAR

Hoof knives with extra fine stainless steel blades.

- Hardwood beech handles with extra fine tip securely riveted to the blade.
- Ideal for exploring holes in the hoof without creating a major excavation.

Sold individually.



Code	Type
220140	Right-handed small blade
220142	Left-handed small blade

FARKNIFE

Carbon steel hoof knives.

- Hardwood beech handles. Long and thick.
- Carbon steel with anti-corrosion treatment.
- More ergonomic than standard hoof knives: grip is particularly suited for farrier working angles.

Sold individually.



Code	Size	Type
220108	L	Right-handed

FARLAND

Carbon steel hoof knives.

- Hardwood beech handles. Long and tapered.
- Carbon steel with anti-corrosion treatment.

Sold individually.



Code	Size
220100	XL

FARLOOP

Carbon steel hoof knives, eye-loop.

- Ambidextrous handle made of polished stained beech.
- The design is ideal for drawing work, with your thumb on the flat end.
- Blades made of carbon steel.

Sold individually.



Code	Size
220110	L
220112	S

HOOFSCRAP

Sharpened hoof scraper with no grip made of shaped steel. Standard quality. Length 30 cm.

Supplied in a plastic case.



Code	Model
220200	Hoofscrap



VETLAND

Hoof knives with carbon steel blades.

- The blade is non-rusting carbon steel and extends through the full length of the composite plastic handle.
- The quality blade sharpens well and retains a good edge.

Sold individually.

HIGHLY
RESISTANT



HOOF TESTERS

Stainless steel hoof testing forceps.

Sold individually.



Code	Model	Type
220002	Single edge	Right-handed narrow blade
220004	Single edge	Right-handed large blade
220006	Single edge	Left-handed narrow blade
220008	Single edge	Left-handed large blade
220010	Double edge (1 short and 1 long)	Right-handed large blade
220032	Double edge (2 long)	Right-handed large blade

Code	Shape	Total length
227100	Flat	350 mm

SHARPENING SET

Sharpening kit for hoof knives and hoof scrapers.

- The kit is composed of 2 items:
 - 6 mm diameter cone-shaped stone for use in an electric drill: ideal for sharpening inside the hook.
 - 10 cm long half-round medium-grit stone: ideal for all the long edges of a hoof knife.
- Supplied in plastic tube.

Sold individually.



Code	Description
220220	Complete set





EQUIPLAST®

Thick layer of highly absorbent cotton wool enclosed in an absorbent cotton gauze cover. High density padding used as a wound dressing to protect against bumps. Great for support, padding, protection and insulation.

- Highly absorbent, it retains any fluid exudate into the dressing.
- The fibres don't stick to the wound as there is a low adherent surface, preserving scar tissues.
- Ideal for equine first aid, stall wraps or shipping wraps.
- Supports and protects fetlocks, hocks and tendons.
- Promotes blood circulation for improved healing.
- Sterilisable.
- Can be cut to the required dimensions.



HIGHLY
ABSORBENT



Code	Size
307280	30 cm x 3.5 m
307282	50 cm x 2 m



A

Adson Brown forceps	40
Adson forceps	40
Allis forceps	41
Animals, large / Equine	69-80
Animals, small	7-36
Anioplast Hydro	21
Anioplast Spongi	21
Anivasc	33
Aniven	32
Apparel	63-64

B

Bags	66
Ball	57
Bandages, cohesive	12
Bandages, elastic	12
Bio Slow feeder bowl	60-61
Blood pressure	47-49
Bottle teats	20
Bouncing ball	57
Bowl, slow feeder	60-61
Bowl, travel	59
Buhner's needle	74
Buhner's vaginal tape	74

C

Caesarean forceps	74
Caesarean scalpel, C-cut	74
Canis Mobile	25
Canis, Handy	25
Carbon tech	44
Carrier, travel	65
Cat, Urolene, urinary catheter	36
Catheters, Tom cat	73
C-cut, caesarean scalpel	74
Cohesive bandages	12
Cooling vest	63
Cosypad, mat	14-15
Cosypop	62
Crazy Urchin Toy	57
Crile forceps	40
Crusher, cutter, tablets	20
Cutter, crusher, tablets	20

D

Danish teat splitter	75
Debakey forceps	40
Dental Toy	56
Digifix, thermometer	18
Digiflash, thermometer	18
Dispoject, safety insulin syringe	34
Dog poop bags	66
Dog poop bags dispenser	66
Dog, Urolene, urinary catheter	36
Doppler, vet	47
Drapvet, surgical drapes	24
Drenching nozzle	77
Drenching syringe	76
Drenching washers	77
Dressing range, Anioplast	21
Drypop	67
Drysoft, puppy pads	19

E

Ecogan, long shoulder gloves	70
E-collars	8-11
E-collars hanger	9
Ecomag	75
Elastic bandages	12
Electrosurgery	50-53
Elplex, hypodermic syringe	72

Enrichment toys, Vétopop	56-58
Equine / Large animals	69-80
Equiplast	81
Esoflex	77
Examed, examination gloves	26
Examination gloves	26
Explogan, long shoulder gloves	70
Extension sets	29

F

Farknife, hoof knives	79
Farland, hoof knives	79
Farloop, hoof knives	79
Feeding syringe	76
Feeding, Vétopop	59-61
Filamid, suture	22
Flat splints	13
Forceps	40-41
Forceps, caesarean	74
Fun Strap, cohesive bandages	12
Funcol	9

G

Garflex, tourniquet	18
Genia gutter splints	13
Geniaplex, hypodermic syringe	72
Gillies needle holder	39
Gloves, examination	26
Gloves, surgical	26
Glycofil, suture	22
Greencol	8

H

Halstead Mosquito forceps	41
Handle, saw wire	75
Handy Canis	25
Hanger, dog E-Collars	9
Hartmann forceps - Alligator	41
Heat mat	44-46
Hoof knives	78-79
Hoof testers	80
Hoofscarp	79
Hook, ovariectomy	41
Hospitalisation, large animals	82
Hospitalisation, small animals	8-21
Hydro, Anioplast, dressing range	21
Hygiene, Vétopop	66-67
Hypodermic syringes	72

I

Imoplast	73
Infuject Site	29
Injection port	28-31
Injection, large animals	72-73
Injection, small animals	28-36
Instruments, surgical	37-41
Insulin syringe, safety, Dispoject	34
Iridectomy scissors	38
IV sets	28-31

K

Kart	25
Killian speculum	39
Knives, hoof	78-79
K-star	29

L

Large animals / Equine	69-80
Lick pad	58
Livestock, large animals	75-79
Long shoulder gloves	70-71

M

Mano Medical	42-49
--------------	-------

Mat, Cosypad	14-15
Mat, Vacuumat, Positioning	16-17
Mathieu needle holder	39
Maxogan, long shoulder gloves	71
Mayo Hegar needle holder	39
Mayo serrated scissors	38
Mayo stille scissors	38
Medical, mano	42-49
Metzenbaum scissors	38
Mobile, Canis	25
Monofil, suture	22

N

Nasoesophageal tubes	77
Nasoflex pump	77
Nasoflex, nasoesophageal tubes	77
Needle holders	39
Needle, rumen	73
Needles, Buhner	74
Needles, Primaject	35
Needles, suture	23
Neocol	9
Nozzle, drenching	77
Nurse watch	18

O

Olsen Hegar needle holder	39
Oralplex	76
Orthopedics, small animals	25
Ovariectomy hook	41

P

Pad, licking	58
Pads, puppy, Drysoft	19
Penmark	76
Perfuvet Duo	28
Perfuvet Micro	29
Perfuvet Plug	28
Perfuvet Twist	30-31
Pet piller	20
Pill popper	20
Pistolman	76
Polyboot	71
Polygon, long shoulder gloves	70
Poop bags	66
Poop bags dispenser	66
Pop lick pad	58
Positioning mat, Vacuumat	16-17
Pressure, blood	47-49
Primaject, needles	35
Progan, long shoulder gloves	70
Protection, large animals	70-71
Protection, small animals	26
puffer, reversible, dog	64
Pump, Nasoflex	77
Puppy pads, Drysoft	19
Puzzle	56

R

Rest, Vétopop	62
Reversible dog puffer	64
Rochester Pean forceps	41
Rumen needle	73
Rumen trocar	75

S

Safety insulin syringe, Dispoject	34
Saw wire	75
Saw wire handle	75
Scissors	38
Screw trocar	75
Seal & Cut	50
Sensigan, long shoulder gloves	71

Set, sharpening	80
Sharpening set	80
Silicone oil, Siligliss	72
Siligliss, silicone oil	72
Slitter, danish teat	75
Slow feeder bowl	60-61
Small animals	7-36
Soft splint roll	13
Spencer scissors	38
Splints	13
Spongi, Anioplast, dressing range	21
Stainless steel, uterine catheter	73
Staple remover	24
Stapvet	24
Stomaflex, nasoesophageal tubes	78
Strap Band, cohesive bandages	12
Surgery, large animals	74
Surgery, small animals	22-24
Surgery, surgical instruments	38-41
Surgical drapes, Drapvet	24
Surgical gloves	27
Surgical instruments	37-41
Surgifar, hoof knives	79
Surgitek, surgical gloves	27
Suture	22-23
Suture needles	23
Suture-rack	23
Syringe, drenching	76
Syringe, feeding	76
Syrman	76

T

Tablet cutter and crusher	20
Teats, bottle	20
Thermometers	18
Tick remover, Tiker	19
Tiker, tick remover	19
Tom cat catheters	73
Tourniquet, Garflex	18
Towel clamp	39
Toy, Crazy Urchin	57
Toy, Dental	56
Training pad, Drypop	67
Transcol	9
Travel bowl	59
Travel pop	65
Treat dispenser puzzle	56
Trocar, rumen	75
Trocar, screw	75

U

Urchin, Toy, Crazy	57
Urinary catheter, cat, Urolene	36
Urinary catheter, dog, Urolene	36
Urolene Cat, urinary catheter	36
Urolene Dog, urinary catheter	36
Uterine catheter, large animals	73

V

Vacuumat, positioning mat	16-17
Vaginal tape, Buhner	74
Velfast	8
Vest, cooling	63
Vet doppler	47
Vetland, hoof knives	80
Vétopop	55-67
Vetrite, examination gloves	26

W

Walking accessories, Vétopop	65
Washers, drenching	78
Watch, nurse	18
Wire, saw	75

NOTES

This image shows a full page of blank, lined paper. It features approximately 20 evenly spaced horizontal blue lines across its entire width. The lines are thin and consistent in color, set against a plain white background. There are no margins, text, or other markings present on the page.

NOTES

[illegible]

HOW TO ORDER YOUR GÉNIA PRODUCTS ?
DIRECTLY THROUGH YOUR WHOLESALER(S)



@GeniaAustralia



Genia Australia



Genia Australia



contact@geniaaustralia.com



www.geniaaustralia.com

génia AUSTRALIA

